

zippie®

Voyage™

**Go Together.  
Grow Together!**



The universal Zippie® Voyage™ early intervention device combines Baby Jogger™ stroller expertise, Zippie's incredibly versatile seating for growing babies, and one of the largest selection of functional accessories in the industry!

**Baby Jogger City Select Base** Zippie teamed up with Baby Jogger to give you an intuitive, easy-to-use stroller option that makes no compromises when it comes to what's best for baby.

**Telescoping Stroller Handle**

Responsive, one-hand steering that adjusts to caregiver heights.

**Canopy with Mesh Window**

When the seat is facing forward, open the mesh window to keep an eye on your little one.



Cross Over



Center Stripe



**Patented Quick-Fold Technology**

Folds into a compact package in one simple step so you can be on your way.



**Puncture-Resistant Rear Wheels**

Travel confidently knowing the Voyage is up for the trip.

**Swivel Front Wheels with Suspension and Lock-Out**

Glide over and around obstacles while keeping baby calm, secure, and comfortable. Lock-out the front wheels for smooth tracking during long-distance strolls.



**Quick Access Wheel Locks**

When juggling kids, bags, or medical equipment, locking the wheels is as simple as flipping one lever on the stroller handle.

**WC-19 Transit Approved**

Approved for use as an occupant's seat in a motor vehicle, the Voyage provides a snug ride for your child.

**Carry-All Storage Basket**

Conveniently store up to 15 lbs. of your essentials.

**Functional Accessories**

The Voyage is available with options and accessories designed to make daily tasks a breeze, so you can delight in the little joys of caring for your family.



IV Support Pole

Cup Holder

Front Mount Utility Platform

Utility Hook

Upper Extremity Support Tray



Folding Booster Base


O2 Cylinder Holder

Battery and Vent Utility Under Frame Tray

Dovetail Vertical Mount (LTV/Trilogy)

Rear Add On Accessory Platform

**Seating for Babies on the Grow!** When it comes to seating, the Voyage team's number one priority is positioning your child for healthy growth and development from babyhood through their toddler years.

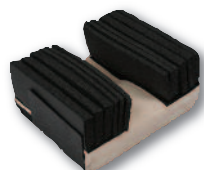

**Adjustable Contour Seat**

- Adjustable Pre-Ischial Ledge
- Lateral Hip Supports
- 3" Seat Depth Growth Kit



**T Back**

- Short (7"– 13")
- Tall (10"–17")
- 1" or 2.5" thick foam




With Hip Blocks Without Hip Blocks

**Standard Planar Seat**

- Contoured foam
- Fixed Ischial Ledge
- Built-in Depth Growth

**Whitmyer Head Supports**

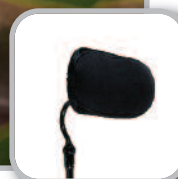

Plush



Adjust-a-Plush



Contour Cradle

**Positioning Accessories**

**Height and Angle Adjustable Footrest** with optional foot straps and shoe holders



Lateral Trunk Support



Butterfly Harness



Medial Support



Foot Straps



Positioning Belt

**Reversible Seat**

Face the seat toward you to interact with your child, or forward so your child can explore their surroundings.



**Tilt/Recline**

Helps improve feeding, digestion, respiratory function, and visual orientation.

**Forward (Anterior) Tilt**

Helps babies achieve developmental milestones such as bringing hands to the mouth, reaching for objects, or engaging in activities.



**Easy Access Tilt/Recline Release**

Tilting and reclining the seat is as simple as engaging bicycle brake-style levers behind the backrest.

# Color Options



**Cross Over Style**

**Center Stripe Style**

## Specifications

### TECHNICAL SPECIFICATIONS

Frame Style:	Early Intervention Device
HCPCS Code:	E1232
	E1234 (without seating)
Chair Weight:	Base only; 23 lbs.
Frame Material:	Aluminum
User Weight Capacity:	60 lbs.
Transit Approved:	Yes (WC-19)

### FRAME AND FEATURE MEASUREMENTS (Dependent on Configuration)

Overall Width:	24"
Seat Widths:	6" – 13"
Seat Depths:	5" – 13"
Back Height:	7" – 17"
Folded Height:	12"
Tilt Range:	Tilt: 50° (+5 to -45°)

### OPTIONS & ACCESSORIES

Frame Type:	Tilt Tilt & Recline
Caster Wheels:	7"
Rear Wheels:	12" Pneumatic with Airless Inserts
Backrest Type:	T-Back Cushion 1" Thick T-Back Cushion 2.5" Thick
Seating:	Standard Planar Seat without Hip Blocks Standard Planar Seat with Hip Blocks Contoured Wedge Seat Contoured Anti-Thrust Seat Adjustable Contour Seat

The Zippie Voyage



How to Unfold & Fold the Zippie Voyage



How to Tilt, Recline, & Reverse the Zippie Voyage Seat



### QR Codes to Voyage Videos