

## Solite Pro Homecare Bed Package

HCPCS:



FIGURE A



FIGURE B



FIGURE C



FIGURE D



FIGURE E

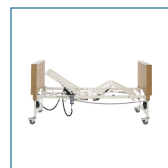


FIGURE F



FIGURE G

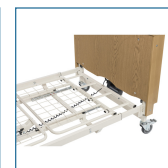


FIGURE H

The Solite Pro from Drive DeVilbiss Healthcare is a fully electrical profiled homecare bed that offers smooth, quiet operation and versatile positioning for maximum patient comfort. The Solite Pro can nest within itself and comes standard with a built-in transport stand for storage and easy transport. Additional features include a tool-free built-in leg extension extending from 80" to 84", knee break function, and improved actuators in the variable height mechanism that add to the strength of the bed, all while retaining an attractive, domestic look and feel.

- Redesigned platform for improved patient comfort when profiling
- Individual sections strengthened for extended performance and reliability
- Improved integral leg extension with safety features to prevent entrapment
- Improved auto-regressing backrest prevents pinching in the sacral area
- All individual sections of the bed weigh less than 55 lbs. (25 kg) for easier assembly
- Four braked castors for maneuverability and safety
- Easy-use handset controls all functions and includes lockout function
- Transport stand included
- Optional integral full-length side rails for safety and ease of use (item #: 15006LEXRAIL)
- Additional optional accessories sold separately include grab handles in two height options, standard (item #: 15006LPSTDRAIL) and short (item #: 15006LSHRTRAIL)
- Rail accessories meet FDA and ASTM entrapment guidelines and standards

PRODUCTS	Description
15006LP	Solite Pro Low Homecare Bed
SPECIFICATIONS	Description
Product Weight	182.1 lbs. (on transport stand)
Bed Height Range	15.5" (mattress platform height) to 32" (mattress platform height)
Casters	3.93" DIA Castors with directional and movement brake (bed must always be used with casters)
Ground Clearance	3.9"
Motor	Linak LA27 Actuators × 4 per bed Stroke lengths HI/LOW: 16.5" stroke, 23.6" install length Head Section: 7.87" stroke, 14.5" install length Foot Section: 2.5" stroke, 10.6" install length
Bed Height Minimum With Casters	11.6"
Bed Height Maximum With Casters	28.2"
ACCESSORIES	Description
15006LPSTDRAIL	Elland standard grab handle, PR
15006LSHRTRAIL	Elland short height grab handle, PR
15006LEXRAIL	Grange telescopic side rails, PR