




ADJUSTER X

*Shown with GlideWear®

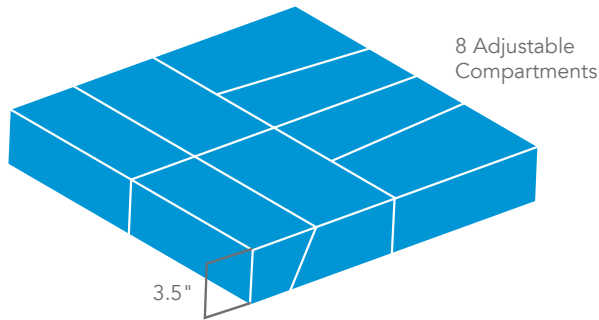

STANDARD SIZES

	16" Wide	18" Wide	20" Wide	22" Wide	24" Wide
16" Deep	16 x 16	18 x 16	20 x 16	22 x 16	24 x 16
18" Deep	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
20" Deep	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20

WEIGHT CAPACITY

All Sizes	550 lbs
-----------	---------


PART NUMBERS

Cover Style	Stretch-Air	Comfort-Tek
Standard	AJX-S- WWDD	AJX-F- WWDD
With GlideWear®	AJX-SQ- WWDD	AJX-FQ- WWDD


WWDD TO BE REPLACED. WW = WIDTH & DD = DEPTH IN INCHES

COVER OPTIONS

STRETCH-AIR

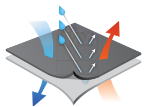

COMFORT-TEK

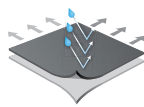

*Available with GlideWear® Technology. See page 3 for details.



Designed with eight adjustable chambers, the Adjuster X is the most comprehensive adjustable skin protection cushion on the market.

Liner is made out of a breathable mesh material and allows the full liner with air cells to be washed and dried. Hot or cold settings may be used.

I.T. bony prominences each have a dedicated chamber that can be partially or fully off-loaded giving the maximum healing opportunity for any stage wound.

Independently adjustable lateral chambers in the rear accommodate lateral pelvic asymmetries or loading the greater trochanters.

Front four chambers can be adjusted to handle seat surface height at front edge variances increasing heel contact for foot propellers, and altering the anti-thrust characteristics for better pelvic positioning.

ADJUSTER



FEATURES

Five fully adjustable compartments provide a custom fit with maintenance free Comfort Cells.

PART NUMBERS

Cover Style	Thickness	Stretch-Air	Comfort-Tek
Standard	3.5"	AJ-S- WWDD	AJ-F- WWDD
Low	2.5"	AJ-SL- WWDD	AJ-FL- WWDD

WWDD TO BE REPLACED. WW = WIDTH & DD = DEPTH IN INCHES

Please refer to the standard size chart in the Adjuster X section for available sizes. Adjuster Low not available in 22" or 24" wide sizes.