

HAULWAY®



Scooter and Power Chair Lifts

HAULWAY lifts allow you to transport your scooter or power chair on the rear of your vehicle safely and efficiently, whether it be across town or on the highway. With easy-to-use features, the HAULWAY takes the worry out of transporting your mobility device with the turn of a key. Each model is simple to install and features a durable and dependable design.



HAULWAY FEATURES

- Key switch activation and easy operation
- Automatic tilt feature keeps your mobility device safe during transport, even on steep inclines
- Folds up for storage when not transporting mobility device
- Models available to carry three- and four-wheeled scooters and most power chairs
- Can be used with class II or III hitch
- Adjustable to allow a working distance of 10-19" from bottom of hitch to ground, which can be increased to 21.5" with available options
- Attractive powder coated finish
- Simple wiring kit and illuminated license plate holder included

Visit ezaccess.com for more information or to locate a dealer near you.

The HAULWAY® Lift Matcher will help you find the lift that best matches your vehicle and mobility device. Visit our website to use this great tool!



HAULWAY Lite

The HAULWAY Lite is a lightweight and easy-to-use lift, designed for smaller three- or four-wheeled scooters, the most popular mobility device category.

A hold down-arm and one retractable ratchet tie-down keep the scooter securely in place. The trays are channeled to hold all wheels of the scooter in position.

Scooter Tray Models Include:

- Lite 3-wheel features a center tray (HLTE3)
- Lite 4-wheel features 2 side trays (HLTE4)
- Lite Deluxe features 3 trays (HLTED)

HAULWAY Standard

The HAULWAY Standard, the most versatile lift in the line, accommodates midsize three- and four-wheeled scooters and power chairs.

A hold-down arm and one retractable ratchet tie-down firmly secure the mobility device to the lift. The scooter trays are channeled to hold all wheels in position and power chair trays are designed for “drive on, drive off” functionality. Add-on trays increase the lift’s versatility and can stand the test of time by providing options that change as needed.

Scooter Tray Models Include:

- Standard 3-wheel features a center tray (HSTD3)
- Standard 4-wheel features 2 side trays (HSTD4)
- Standard Deluxe features 3 trays (HSTDD)

Power Chair Tray Models Include:

- Standard Power Chair features 2 side trays (HSTDPC)
- Standard Power Chair Deluxe features 3 trays (HSTDPCD)



HAULWAY Heavy Duty

The HAULWAY HD is a heavy-duty lift featuring the highest weight capacity in the line, making it ideal for midsize to large scooters and power chairs.

In place of a hold-down arm, the lift features a durable and easy-to-use 4-point retractable ratchet system to ensure secure transport. The tray is designed for “drive on, drive off” functionality.

Scooter and Power Chair Tray Models Include:

- Heavy Duty Universal features 1 large tray (HHD)

Quick Specs

	Weight	Weight Capacity	Width x Length	Hitch	Ideal for:
Lite 3-wheel	85 lbs.	175 lbs.	22" x 43"	II or III*	Compact 3-wheeled scooter
Lite 4-wheel	89 lbs.	175 lbs.	22" x 43"	II or III*	Compact 4-wheeled scooter
Lite Deluxe	93 lbs.	175 lbs.	22" x 43"	II or III*	Compact 3- or 4- wheeled scooter
Standard 3-wheel	90 lbs.	275 lbs.	26.5" x 45"	II or III*	3-wheeled scooter
Standard 4-wheel	98 lbs.	275 lbs.	26.5" x 45"	II or III*	4-wheeled scooter
Standard Deluxe	106 lbs.	275 lbs.	26.5" x 45"	II or III*	3- or 4- wheeled scooter
Standard Power Chair	107 lbs.	275 lbs.	26.75" x 38"	II or III*	Power chair or 4-wheeled scooter
Standard Power Chair Deluxe	112 lbs.	275 lbs.	26.75" x 45"	II or III*	Power chair or 3- or 4-wheeled scooter
Heavy Duty Universal	107 lbs.	325 lbs.	27" x 45"	II or III*	Medium to large power chair or scooter


Lite



Standard



Heavy Duty

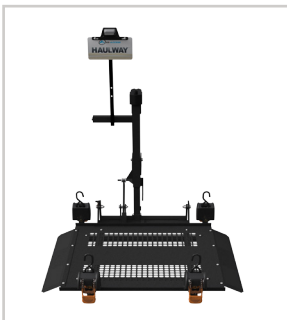

* Please consult your HAULWAY dealer to determine your ideal hitch and lift combination.

OPTIONS



Swingaway



Battery Pack



Scooter and Power Chair Cover

- The HAULWAY Swingaway enables full access to the vehicle's tailgate, rear door, or cargo area, as it allows the lift (without mobility device) to be swung out of the way. The Swingaway option can be used with all HAULWAY models. Available with a class III hitch.
- The HAULWAY portable battery pack allows for increased ease of installation and a Plug-N-Play solution. Wiring directly to the vehicle's battery is not required when this convenient option is chosen. The battery pack is simply installed near the license plate holder on a provided bolt-on tray.
- The Scooter and Power Chair Covers from our EZ-ACCESSORIES® line are water-resistant and provide protection from dust as well as the sun. Visit our website to view additional accessories designed for power mobility devices.



Class II No Drop



Class II 2" Drop



Class II 4" Drop



Class III No Drop



Class III 2" Drop



Class III 4" Drop

- Class II and III hitch adapters (L-tubes) are available up to a 4" drop. The type of L-tube or Swingaway needed is dependent on the distance from the bottom of the vehicle's hitch receiver to the ground.
 - » If the distance is 10"-17.5", a standard (no drop) L-tube or Swingaway is required.
 - » If the distance is 17.5"-19.5", a 2" drop L-tube or Swingaway is required.
 - » If the distance is 19.5"-21.5", a 4" drop L-tube or Swingaway is required.

All HAULWAY Scooter and Power Chair Lifts are covered under the EZ-ACCESS 3-year warranty and are proudly made in the USA.



customerservice@ezaccess.com | 1.800.451.1903 | www.ezaccess.com