

**Jazzy®**  
600/600XL

TRUE ALL-AROUND PERFORMANCE.  
MAXIMUM STABILITY. NO COMPROMISES.

## Jazzy® 600/600XL

The Jazzy® 600 and 600XL featuring Mid-Wheel 6® utilize Pride's original, patented mid-wheel drive with enhanced stability. Mid-Wheel 6 incorporates ATX Suspension for extra stability with independent active front OMNI-Casters (nylon, spherical-shaped casters) that raise and lower as needed to compensate instantly for changes in the terrain. They prevent hang-ups and allow for stable transitions from one surface to the next, regardless of the angle of approach.

Mid-Wheel 6® provides superb all-around performance with the added benefit of 6-wheel stability. Whether transitioning curbs, climbing rough terrain, or reverse handling, Mid-Wheel 6 consistently excels.



**Pride®**  
Mobility Products Corp.

800-800-8586 (US) • 888-570-1113 (Canada) • [www.pridemobility.com](http://www.pridemobility.com)



The Jazzy® 600 series models include:

- Jazzy 600 2S-C
- Jazzy 600 2S-SS
- Jazzy 600 2SP-C (250 lbs. weight capacity)
- Jazzy 600 2SP-SS (250 lbs. weight capacity)
- Jazzy 600 3S-C
- Jazzy 600 3S-SS
- Jazzy 600 3SP-C (250 lbs. weight capacity)
- Jazzy 600 3SP-SS (250 lbs. weight capacity)
- Jazzy 600 XL-C (450 lbs. weight capacity)
- Jazzy 600 XL-SS (450 lbs. weight capacity)

## FEATURES

- OMNI-Casters (nylon, spherical-shaped casters) on front and rear to prevent wheel hang-ups
- Side-mounted, easily accessible freewheel levers
- Active-Trac® ATX Suspension incorporates front OMNI-Casters for enhanced performance over more varied terrain
- Easy front access to batteries
- Easily accessible side charger port


## OPTIONS

- Solid seat pan
  - Cup holder
  - Swing-away joystick
  - Swing-away footrests
  - Elevating leg rests
  - Oxygen tank holder
  - Cane/crutch holder
  - Walker holder
  - TRU-Balance® Power Tilt (N/A on XL)
  - Synergy® Seat (400 lbs. maximum weight capacity)
- <sup>1</sup> 22" x 20" Contoured seat  
20" x 20" with Synergy® Seat or TB<sup>2</sup> Tilt

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3 %.

	Jazzy® 600	Jazzy® 600XL
Weight capacity:	300 lbs.	450 lbs.
Maximum speed:	Up to 5 mph	Up to 4 mph
Ground clearance:	3"	3"
Turning radius:	20"	20"
Overall length:	36.5" (w/o front riggings)	37.5" (w/o front riggings)
Overall width:	23"	23"
Wheels:	Drive: 14"	Drive: 14"
Front:	OMNI-Casters: 5" solid	OMNI-Casters: 5" solid
Rear:	OMNI-Casters: 6" solid	OMNI-Casters: 6" solid
Suspension:	Active-Trac® ATX Suspension	Active-Trac® ATX Suspension
Drive train:	Two-motor, Mid-Wheel 6°	Two-motor, Mid-Wheel 6°
Braking system:	Intelligent braking (electronic, regenerative disc brakes)	Intelligent braking (electronic, regenerative disc brakes)
Std. Electronics:	60 amp, PG VR2 controller	70 amp, PG VR2 controller
Specialty controls:	Yes	N/A
Battery charger:	On-board, 4 amp (standard)	On-board, 4 amp (standard)
Per-charge range:	Up to 20 miles	Up to 15 miles
Battery requirements:	(2) 12 volt, deep cycle	(2) 12 volt, deep cycle
Battery size:	NF-22	NF-22
Transport weight		
Battery weight:	38 lbs. each	38 lbs. each
Base weight:	127 lbs.	133.5 lbs.
Std. seat weight:	37 lbs. (medium-back)	42.5 lbs. (high-back)
Max. Seat dimensions:	22" W x 20" D <sup>1</sup>	24" W x 24" D
Warranty:	Lifetime limited warranty on frame; 2-year limited warranty on electronics; 18-months limited warranty on drive motors	

Colors:

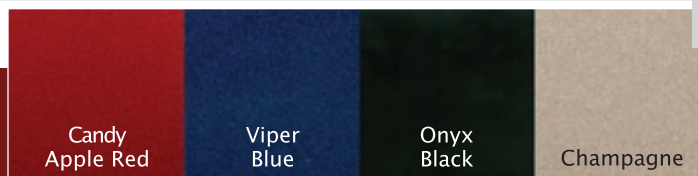

Candy  
Apple Red


Viper  
Blue

Onyx  
Black

Champagne

## Pride SilverStar Lifts & Ramps

- Backpacker® Plus
- Full Platform
- Power chair lift with straps

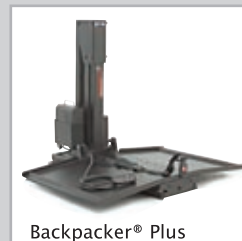

Backpacker® Plus



Full  
Platform

These lifts are compatible with the Jazzy® 600/600XL.

800-800-8586 (US) • 888-570-1113 (Canada) • www.pridemobility.com

**Pride**<sup>®</sup>  
Mobility Products Corp.