



ROHO[®] SOFFLEX[®] 2 MATTRESS OVERLAY SYSTEM

SKIN PROTECTION

Protecting an individual's skin/soft tissue is critical for prevention and treatment of pressure ulcers. The ROHO SOFFLEX 2 Mattress Overlay System is a non-powered, adjustable and reactive support surface for individuals who require skin/soft tissue protection or have other therapeutic needs. The SOFFLEX 2 provides individualized pressure, shear and friction management through force redistribution to an individual's shape.

COST EFFECTIVE

The SOFFLEX 2 is a cost effective option for different care settings, costing less than high-quality powered mattress systems. It can also be easily cleaned and disinfected for reuse.

ADJUSTABLE

Flexible, interconnected air cells can be adjusted to match the body's contour, providing personalized therapy and critical skin/soft tissue protection. The SOFFLEX 2 is made of three sections that snap together. Simply add or release air in each mattress overlay section to accommodate the individual's needs.

NON-POWERED

No need for concern of power failure or energy consumption costs.

EASE OF USE

The SOFFLEX 2 is lightweight, easy to transport, easy to adjust and set up on different bed types.

Choose the SOFFLEX 2 for protection and peace of mind!



ROHO® SOFFLEX® 2 MATTRESS OVERLAY SYSTEM

APPLICATIONS

THE SOFFLEX 2 IS EFFECTIVE FOR INDIVIDUALS WHO:

- are at risk of developing pressure ulcers and/or tissue injury
- are at risk for skin/soft tissue breakdown
- have any stage pressure ulcer, skin grafts or recent myocutaneous flap surgery
- need comfort
- need pain management
- need a mattress overlay that provides an individualized fit

SPECIFICATIONS

INCLUDES:

SOFFLEX 2 Mattress Overlay System (3 sections), hand inflation pump, straps, operating instructions and repair kit.

CONSTRUCTION:

Flame-retardant Polyurethane

SIZES:

MATTRESS OVERLAY SYSTEM (3 Sections)

SYSTEM: 36" x 81.75" x 3.5" (91.5 cm x 207.5 cm x 9 cm)

SECTION: 36" x 27.5" x 3.5" (91.5 cm x 70 cm x 9 cm)

WEIGHT:

MATTRESS SYSTEM: 6.5 lbs (2.9 kg)

SECTION: 2.25 lbs (1 kg)

WEIGHT LIMIT:

330 lbs (150 kg) for complete system

WARRANTY:

18-month limited warranty

U.S. MEDICARE PART B CODE:

E0371 Approved for Group II Support Surfaces

FEATURES AND BENEFITS

AIR FLOATATION TECHNOLOGY:

Air Floatation technology redistributes forces that can cause peak pressure, shear and friction. The SOFFLEX 2 is effective for comfort, pain management, skin/soft tissue protection and treatment.

TRANSPORTABLE:

The SOFFLEX 2 weighs only 6.5 lbs (2.9 kg). The air can be released for easy transport between the hospital, care facility, home or travel. No need to buy a new mattress overlay for each setting.

ADJUSTABLE:

Three sections create a multi-zoned support surface for an individualized fit. Each mattress overlay section can be adjusted per body segment to address the specific needs of the individual's shape and size.

FITS MOST BEDS:

Can be placed on top of almost any size mattress, including the individual's existing home mattress, hospital bed mattress or care facility bed mattress.

COST-EFFECTIVE:

The SOFFLEX 2 is non-powered, reusable and costs less than powered mattress systems.

LIGHTWEIGHT:

A complete SOFFLEX 2 system weighs only 6.5 lbs (2.9 kg). Easy to carry and handle during set up, inflation and adjustment.

NON-POWERED:

Energy-efficient because it requires no power source for use. There is also no noise or heat from a pump.

COVER OPTIONS:

Two options designed to fit every ROHO Mattress Overlay System: ROHO® Reusable Mattress Overlay Cover (one piece that tucks in like a fitted sheet) or the ROHO® Mattress Overlay Enclosure Cover (has a top and bottom that zip together to completely enclose the mattress overlay). Each cover option is fluid-resistant, flame-retardant, and can be easily removed for cleaning and disinfecting between uses.

MORE PRODUCTS FROM ROHO



ROHO® Hybrid Elite® Cushion




ROHO® QUADTRO SELECT® HIGH PROFILE® Cushion



ROHO® MOSAIC® Cushion



ROHO® HEAL PAD® Cushion

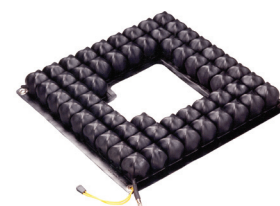

ROHO® Shower/Commode Cushion

FOR ADDITIONAL INFORMATION, PLEASE CONTACT:



100 North Florida Avenue • Belleville, IL 62221-5429 • USA • www.therohogroup.com

U.S.: 1 800 851 3449 • fax 1 888 551 3449

Outside the U.S.: 1 618 277 9150 • fax 1 618 277 6518

©2002, 2010, 2011 ROHO, Inc. The following are registered trademarks of ROHO, Inc.:

ROHO®, HEAL PAD®, Hybrid Elite®, HIGH PROFILE®, MOSAIC®, QUADTRO SELECT®, Shape Fitting Technology® and SOFFLEX®.

ROHO, Inc. has a policy of continual product improvement and reserves the right to amend this document. The current version of this document is available on our website at www.therohogroup.com.

LITD00258 Rev. 05/11 05/11

