

POWER CHAIRS



QUICKIE®


Contents

01	S-646 SE
02	S-646
03	S-525
04	P-222 SE
05	P-220
06	P-200
07	Freestyle
08	G-424
09	Z-500
10	V-521
11	V-121
12	V-100
13	Optional Controls
14	Seating Options
15	Power Matrix (Images)
16	Power Matrix (Chart)

FULL FEATURE REHAB POWER BASE WITH EXTRA POWER

The Quickie S-646 SE offers innovative performance as well as versatile seating and positioning, four pole motors and a top speed of 8.5 mph. PosiDrive Technology™, FoxShox™ suspension and soft roll casters provide a smooth, intuitive ride. Quickie SmartSeat®, TRAX seating and a CG Tilt option allow the Quickie S-646 SE to accommodate general and specialized rehabilitation needs. Form, fit and function with a punch — the S-646 SE has it all.

S-646 SE

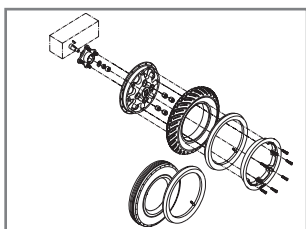

POSIDRIVE TECHNOLOGY

Our new rock solid drive train increases the overall strength of the chair, features easy-to-service hubs and maximizes the chair's overall performance. (Featured on the S-646, S-646 SE, P-222 SE, P-220 and the S-525.)



QUICKIE SMARTSEAT

The Quickie S-646 and S-646 SE can be ordered with a SmartSeat — a breakthrough in tilt and enhanced recline technology. Its unique, patented design offers all the benefits of tilt and recline while preserving posture and skin integrity.



FREEWHEEL LEVER

Easiest free-wheel push force in the industry! The re-designed gear releases are now easier to access and designed to completely disengage the motor, making the chair very maneuverable for transfers.

INTEGRATED POWER C. G. TILT

Integrated sliding center-of-gravity tilt offers up to 53° of tilt. The Quickie S-646 and S-646 SE with tilt can even be driven when tilted to 15°, and it accommodates riders weighing up to 400 lbs.



S-646

FULL FEATURE REHAB POWER BASE

The Quickie S-646 offers all the great features found in the S-646 SE with two pole motors and a top speed of 7mph. PosiDrive Technology™, FoxShox™ suspension and soft roll casters provide a smooth, intuitive ride. Quickie SmartSeat®, TRAX seating and a CG Tilt option allow the Quickie S-646 to accommodate general and specialized rehabilitation needs.

CUSTOMIZED REHAB POWER BASE

The modular design of the Quickie S-525 makes it the ideal chair for those with special seating needs. The rider can choose from any one of sixteen different seat height/angle combinations to get the fit that's right for their individual positioning needs. Multiple performance options and multiple frame sizes further ensure that the chair suits the rider's body and lifestyle.



INTEGRATED POWER C.G. TILT

Integrated, sliding center-of-gravity tilt offers up to 50° of tilt while maintaining a low 18.5" seat-to-floor height. Available to retro fit.



NEW STANDARD WHEEL OPTION

The NEW split rim wheel ships standard on all Quickie S-525 models, which now feature our state of the art Posi-Drive Technology™ plus a new 12"x3" tire for greater terrain hugging traction.



P-222 SE

DYNAMIC POWER BASE WITH EXTRA POWER

The powerful and compact design of the Quickie P-222 SE offers the optimum in performance, style and maneuverability. New four pole motors provide a top speed of 8.5 mph. The sleek, compact styling is visually appealing and allows the rider to maneuver in tight environments. The quick-release power pack and optional manual conversion kit deliver both convenience and transportability.



MANUAL WHEELCHAIR CONVERSION KIT

Manual or Power? It's your choice. With the Quickie P-222 SE and P-220 the rider can switch to a manual chair without having to change the seating.



QUICK-RELEASE POWER PACK

Transportation and service is quick, easy and convenient with the removable power pack.



FREE-WHEEL RELEASE LEVER

The free-wheel release levers for the Quickie P-222 SE and Quickie P-220 have been designed for easy access and located where they are protected from impacts. Once the free-wheel mode has been activated minimal push force is required to maneuver the chair.



INTEGRATED POWER C. G. TILT

Integrated sliding center-of-gravity tilt offers up to 50° of tilt. The Quickie P-220 and P-222 SE with tilt can even be driven when tilted to 15° and accommodate riders weighing up to 250 lbs.



P-220

DYNAMIC POWER BASE

The Quickie P-220 is Quickie's price and performance leader. The sleek design, coupled with engineering features such as PosiDrive Technology™, 6.5 mph top speed in its class and our new integral controller, sets the Quickie P-220 apart in terms of quality and performance. The Quickie P-220 delivers superior function, performance and style to meet a variety of mobility needs.

P-200

CLASSIC PERFORMANCE POWER BASE

The Quickie P-200 is the time-proven standard when it comes to responsive acceleration and “turn-on-a-dime” maneuverability. The modular design of the power drive system also makes it possible for the Quickie P-200 to convert quickly and easily to manual operation.



MANUAL CONVERSION

The easy-to-use manual conversion kit is a great way to maintain seating while utilizing manual means of mobility.



POWER DRIVE SYSTEM

The Quickie P-200 was dubbed by Popular Mechanics as the chair that changed everything. The power drive system is still popular in the market today.



INTEGRATED POWER C.G. TILT
Integrated, sliding center-of-gravity tilt offers up to 50° of tilt while maintaining a low 18.5" seat-to-floor height. Available to retro fit.



CAPTAIN'S SEAT
The Quickie Freestyle offers comfortable options for seating.

FREESTYLE MINI-M11
Also available in Full Size.



Freestyle

IMAGINE THE POSSIBILITIES

The unique independent all-wheel suspension and 3-bar linkage system are integral parts of the Quickie Freestyle "Overland Suspension". The Quickie Freestyle features one of the smallest footprints in the industry for a mid-wheel drive chair. This means that the Freestyle can maneuver in tight confined living quarters and fit into small storage spaces that larger mid-wheel drive chairs are unable to manage.





OPTIONAL INTEGRATED POWER C.G. TILT

Get up to 50° of tilt, as much as 6° of pre-tilt and seat height adjustability with the integrated power tilt option. (Available as a retro fit kit)



EASY BATTERY ACCESS

Batteries are easily accessed while rider is seated in the Quickie G-424.



G-424

CENTER DRIVE POWER BASE

The center-wheel drive design of the Quickie G-424 delivers an extremely maneuverable and stable ride. The tight turning radius and a 25" width allow access to confined areas. Forward pitching and rocking are eliminated with the articulating, constant-contact casters and the unique six-wheel suspension. The color-infused body is extremely durable and will look great for years to come.

PEDIATRIC POWER BASE

The Quickie Z-500 provides the optimum power wheelchair base for pediatric seating and positioning making it ideal for the young rider. The unique seat and base design allows for multiple growth adjustments. Seat growth kits are also available. This model allows young riders to explore environments independently.

Z-500



Option



QUICKIE Z-500 & KID KART® XPRESS SEATING

Can be configured to accept Kid Kart Xpress seating for the easy switch between manual and powered mobility.



QUICKIE Z-500 POWER TILT

Get up to 30° of power tilt with our TK Power Tilt System with through drive control option.

V-521

MODERATE DUTY FOLDING CHAIR

The Quickie V-521 offers a unique combination of style, portability and function. Reinforced, oval aluminum tubing makes the design both attractive and extremely durable. The battery pack can be quickly removed, allowing the chair to be easily folded and stowed for travel. This robust, transportable chair provides a power option for any budget.



QUICK-RELEASE BATTERY CAGE

Transporting the Quickie V-521 is easy and convenient. Simply raise the handle, and the battery case lifts out of the chair.



FOLDED CONFIGURATION

Regardless of the seat width setting, the Quickie V-521 and V-121 fold to less than 13" wide.



CONVENIENT AND EASY FREE-WHEELING

Easy to use and easy to reach. A quick turn of the lever is all it takes to engage or disengage free-wheeling on both the Quickie V-521 and V-121.



QUICK-RELEASE BATTERY CAGE

Disassembly and transport is made easy thanks to our removable battery cage with handle.



V-121

LIGHT DUTY FOLDING CHAIR

The Quickie V-121 is designed for the rider. The simplicity of this model makes it very affordable and easy to operate. The batteries can be quickly removed, and the chair can then be folded for easy transport.

V-100

LIGHT DUTY RIGID CHAIR

The Quickie V-100 combines performance and strength in a lightweight transportable chair. The rigid frame offers a durable, solid ride and a stable structure for seating and positioning accessories. Fold the seat back down for easy transport.



CONVENIENT AND EASY FREE-WHEELING
Easy to use and easy to reach. A quick turn of the lever is all it takes to engage or disengage free-wheeling.



FOLD-DOWN ANGLE-ADJUSTABLE BACKREST
Angle-adjustable backrest accommodates posture and positioning systems. Folds down for maximum portability. (Shown without batteries)

Optional Controls

THE POWER OF CHOICE

Quickie Power Series can be controlled with a variety of specialty input devices. The Quickie QTRONIX® system provides customized control for virtually any rider. Refer to the Quickie Power Feature Matrix at the back of this brochure for more information regarding optional control and chair compatibility. Reference the Quickie Power Feature Matrix for compatibility. For more information on these and other specialty control options, please call customer service.

QTRONIX® SYSTEM



Power Control Module



Universal Specialty Control Module



Remote Joystick - 5 program options



Programming Pad



Actuator Lights Module (ALM)

MULTI-SWITCH INPUTS



Zero Touch Switches (5)



Breath Control



Star Board



Proximity Head Array



Actuator Lights Module Wafer Board(ALM)



Penta Switch



Tri-Switch Head Array

PROPORTIONAL INPUTS



Attendant Control



Chin Control



Proportional Head Control



Mini Proportional Joystick

SINGLE SWITCH INPUTS



Soft Switch



Treadle Switch



Grip Switch



Ribbon Switch



Buddy Button



Micro-Switch

Not Shown:

Zero Touch Switches
(Tray Mount)

CA-5 Interface Box

Wobble Switch

Disc Switch

TM-40 Power Seating
Module

Seating Options

CUSTOMIZED SEATING

The Quickie Power Series can be customized with a choice of seating options to offer the rider the ultimate in comfort and support.



HIGHBACK CONTOUR (G-424, S-525, FREESTYLE)

The Highback Contour seat optimizes support and allows for easy lateral transfers. The contour foam maximizes comfort and virtually eliminates bottoming out.



CONTOUR SEAT (G-424, S-525)

The molded foam construction of the Contour Seat provides both support and comfort.



SEAT FRAME (G-424, S-525, P-222 SE, P-220, FREESTYLE)

The Seat Frame seat is available in a variety of sizes to ensure a perfect fit for riders with more specialized seating needs.



QUICKIE TRAX™ SEAT (S-646, S646 SE AND SMARTSEAT)

The sliding track mounts of the Quickie TRAX™ Seat allow for convenient seat and armrest adjustments and a wide range of seating positions.

S - SEATING AND POSITIONING



SmartSeat: Power Tilt and Enhanced Recline



S-646 SE: Full Feature Rehab Power Base



S-646: Full Feature Rehab Power Base



S-525: Standard Rehab Power Base

P - PERFORMANCE SERIES



Freestyle: Dynamic Power Base



P-222 SE: Dynamic Power Base



P-220: Classic Performance



P-200: Classic Performance

G - GENERAL MOBILITY



G-424: Center Drive Power Base

Z - PEDIATRIC SERIES



Z-500: Moderate Duty Pediatric Power Base

V - VALUE SERIES



V-521: Moderate Duty Folding Chair



V-121: Light Duty Folding Chair



V-100: Light Duty Rigid Chair

Power Matrix



TRANSIT TIE-DOWNS

--	--	--	--	--	--	--	--	--	--	--	--	--

SEAT TYPES

Folding Seat Frame												
Rigid Seat Frame												
Quickie Trax Seat Frame												
Contour												
Highback Contour												
SmartSeat												

SEAT FRAME DIMENSIONS

Width	14"-24"	14"-24"	14"-20"	14"-24"	14"-24"	14"-20"	14"-24"	14"-20"	12"-16"	14"-20"	14"-20"	14"-20"
Depth	14"-22"	14"-22"	14"-20"	14"-20"	14"-20"	14"-20"	14"-20"	14"-20"	10"-15"	14"-18"	14"-18"	14"-18"
Back Height	15"-24"	15"-24"	15.5"-19.5"	15.5"-19.5"	15.5"-19.5"	15.5"-19.5"	15.5", 17.5", 19.5"	15.5"-19.5"	11"-16"	15"-20"	16"-19"	15"-20"
Back Angle (°)	84°-124°	84°-124°	86°-102°	86°-102°	86°-102°	56°-102°		86°-102°		0 (85°-115°)	90°-98°	0 (85°-95°)
Seat-to-Floor Height-Front (non-tilt systems)	18"-21"	18"-21"	18"-21"	18"-21"	18"-21"	18.5"-19.5"	15.5"-21.5"	18"-21"	20"	18.5"/20.25"	18.5"	18.5"/20"
Seat-to-Floor Height-Front (tilt systems)	18.5"/20.5"	18.5"/20.5"	18.5"-20"	19.75"-21.75"	19.75"-21.75"		17.5"-23"	18.5"-20"				

POWER SEATING

Manual Tilt												
Power Tilt												
SmartSeat Power Recline												
SmartSeat Power Tilt & Recline												

WARRANTY

Frame	Lifetime	Lifetime	Lifetime	Lifetime	Lifetime	Lifetime	Lifetime	Lifetime	Lifetime	Lifetime	Lifetime	5 Year
Integral Controller	NA	NA	2-yrs	NA	2-yr	NA	2-yrs	2-yrs	NA	2-yrs	2-yrs	2-yrs
Remote Joystick	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr			
Remote Control Module	2-yrs	2-yrs	2-yrs	2-yrs	2-yr	2-yr	2-yrs	2-yrs	2-yrs			
Motor	2-yrs	2-yrs	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr
Gearbox	2.5-yrs	2.5-yrs	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr	1-yr

* Reduce overall width by 2.5" when using the optional 12" X 3" drive tire.

** Add 2.5" to the overall width when using the optional 12" X 4" drive tire.

O - Option

S - Standard

Power Matrix

SELECT THE CHAIR THAT IS RIGHT FOR YOU



	S-646 SE	S-646	S-525	P-222 SE	P-220	P-200	Freestyle	G-424	Z-500	V-521	V-100	V-121
POWER BASE												
Frame												
MAXIMUM SPEED (MPH)												
Standard	8.5	7	4.5	8.5	6.5	7	5.5	4.5	5	4.5	5	4.5
Option	7	6	6.5				6	6.5				
DRIVE WHEELS												
Standard Size	14 x 3	14 x 3	12 x 3	12 x 4	12 x 4	10 x 3	10 x 3	14 x 3	12 x 2	12 x 2	12 x 2	12 x 2
Optional Size			12 x 4	12 x 3	12 x 3		14 x 3					
CASTER SIZE												
8"	•	•	•	•		•		•	•	•	•	•
9"	•	•		•								
REAR SUSPENSION												
	•	•	0	•			S	•				
FRONT SUSPENSION												
	0	0	0	0	0		S	•				
Overall Width												
Standard	25.75	25.75	22.75**	24.5*	24.5*	24	22.5	25.0	22.75	19.5-25.5	22.5-26.5	19.5-25.5
Option 1	27.75	27.75	24.75**	26.5*	26.5*	26	24.5					
Option 2			26.75**									
CONTROL SYSTEM												
Integral			S		S		S	S		S	S	S
OTRONIX Remote Joystick	S	S	0	S	0	S	0	0	S			
Specialty Controls	0	0	0	0	0	0	0	0	0			
Battery Type												
U1			0				0		0	0	0	0
22NF	0	0	0	0	0		0	0	0	0	0	0
Gp 24	0	0		0	0	0	0	0				
MAX WEIGHT CAPACITY												
150 lbs									S			
250 lbs			S			S		S			S	S
300 lbs					S		S			S		
350 lbs			0	S				0				
400 lbs	S	S										

We highly recommend distribution through a certified Rehab Professional when fitting and ordering any Quickie Products.

QUICKIE®