

\* HCPCS CODING:  
E2G13 / E2G14

Positioning Back




# ACTA-BACK


\*see chart for dimensions

## MOLDABLE CENTER SECTION

- Backrest pad may be removed and molded, increasing support and comfort to accommodate lumbar, scoliosis, and kyphosis needs.

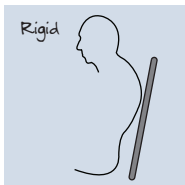


- Moldable aluminum stays are offset from the center, removing pressure points from the spine.

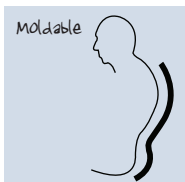

Patented Design: #6,753,074

## - NEW FEATURES -

- Lightweight aluminum shell.
- Provides hip and scapula relief.
- Three layer foam construction, including Visco topper layer.
- Laterals and headrests can be attached to all Acta-Back models (see pages 22 and 23).

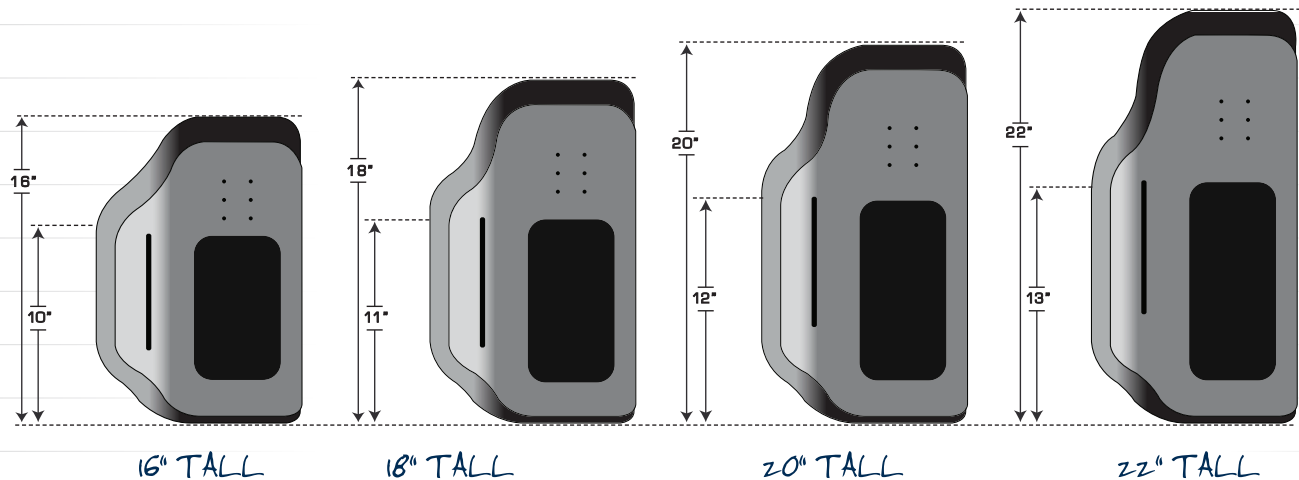



Quick Release Feature

## - COMPASS HARDWARE -

Compass Hardware: Patented #7,104,610

- Patented Compass Hardware is fully adjustable for height, depth, and angle of recline.
- Requires only 1.5" of space on each cane for mounting purposes.
- All adjustments can be made while the user is seated in the wheelchair, eliminating the need for multiple transfers during the fitting process.
- Quick-release feature allows backrest to be removed vertically within seconds.
- Hardware may be installed within minutes and fits 3/4", 7/8", and 1" tubing.
- New positive locking feature for angle rotation.

## - COVERS -

**STRETCH-AIR FABRIC** — designed for air transmission and heat reduction. This breathable laminate combines a high-stretch polyester top layer with a waterproof, breathable vapor barrier.

**COMFORT-TEK** — fluid proof exterior has supple leathery properties that provide an easily cleaned surface and protection for the inner cushion. Comfort-Tek has multidimensional stretch and helps alleviate pressure on bony prominences.

BACK HEIGHT: 10" Tall		FABRIC TYPE: <i>Stretch-Air</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-MLUC-14</b>	9.5"	2.8 lbs.	13 - 15	2.5"	250 lbs.	
<b>BK-MLUC-16</b>	9.5"	2.9 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-MLUC-18</b>	9.5"	3.1 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-MLUX-20</b>	9.5"	3.2 lbs.	19 - 21	2.5"	300 lbs.	

BACK HEIGHT: 16" Tall		FABRIC TYPE: <i>Stretch-Air</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-MUC-16</b>	16"	4.4 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-MUC-18</b>	16"	4.7 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-MUC-20</b>	16"	5.0 lbs.	19 - 21	2.5"	300 lbs.	
<b>BK-MUC-22</b>	16"	5.3 lbs.	21 - 23	2.5"	350 lbs.	
<b>BK-MUX-24</b>	16"	5.8 lbs.	23 - 25	2.5"	400 lbs.	

BACK HEIGHT: 18" Tall		FABRIC TYPE: <i>Stretch-Air</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-M18UC-16</b>	18"	4.8 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-M18UC-18</b>	18"	5.2 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-M18UC-20</b>	18"	5.4 lbs.	19 - 21	2.5"	300 lbs.	
<b>BK-M18UC-22</b>	18"	5.7 lbs.	21 - 23	2.5"	350 lbs.	
<b>BK-M18UX-24</b>	18"	6.3 lbs.	23 - 25	2.5"	400 lbs.	

BACK HEIGHT: 20" Tall		FABRIC TYPE: <i>Stretch-Air</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-M20UC-16</b>	20"	5.2 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-M20UC-18</b>	20"	5.5 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-M20UC-20</b>	20"	5.9 lbs.	19 - 21	2.5"	300 lbs.	
<b>BK-M20UC-22</b>	20"	6.2 lbs.	21 - 23	2.5"	350 lbs.	
<b>BK-M20UX-24</b>	20"	6.8 lbs.	23 - 25	2.5"	400 lbs.	

BACK HEIGHT: 22" Tall		FABRIC TYPE: <i>Stretch-Air</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-MTUC-16</b>	22"	5.6 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-MTUC-18</b>	22"	5.8 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-MTUC-20</b>	22"	6.2 lbs.	19 - 21	2.5"	300 lbs.	
<b>BK-MTUC-22</b>	22"	6.5 lbs.	21 - 23	2.5"	350 lbs.	
<b>BK-MTUX-24</b>	22"	7.0 lbs.	23 - 25	2.5"	400 lbs.	

BACK HEIGHT: 10" Tall		FABRIC TYPE: <i>Comfort-Tek</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-MLFC-14</b>	9.5"	2.8 lbs.	13 - 15	2.5"	250 lbs.	
<b>BK-MLFC-16</b>	9.5"	2.9 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-MLFC-18</b>	9.5"	3.1 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-MLFX-20</b>	9.5"	3.2 lbs.	19 - 21	2.5"	300 lbs.	

BACK HEIGHT: 16" Tall		FABRIC TYPE: <i>Comfort-Tek</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-MFC-16</b>	16"	4.4 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-MFC-18</b>	16"	4.7 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-MFC-20</b>	16"	5.0 lbs.	19 - 21	2.5"	300 lbs.	
<b>BK-MFC-22</b>	16"	5.3 lbs.	21 - 23	2.5"	350 lbs.	
<b>BK-MFX-24</b>	16"	5.8 lbs.	23 - 25	2.5"	400 lbs.	

BACK HEIGHT: 18" Tall		FABRIC TYPE: <i>Comfort-Tek</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-M18FC-16</b>	18"	4.8 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-M18FC-18</b>	18"	5.2 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-M18FC-20</b>	18"	5.4 lbs.	19 - 21	2.5"	300 lbs.	
<b>BK-M18FC-22</b>	18"	5.7 lbs.	21 - 23	2.5"	350 lbs.	
<b>BK-M18FX-24</b>	18"	6.3 lbs.	23 - 25	2.5"	400 lbs.	

BACK HEIGHT: 20" Tall		FABRIC TYPE: <i>Comfort-Tek</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-M20FC-16</b>	20"	5.2 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-M20FC-18</b>	20"	5.5 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-M20FC-20</b>	20"	5.9 lbs.	19 - 21	2.5"	300 lbs.	
<b>BK-M20FC-22</b>	20"	6.2 lbs.	21 - 23	2.5"	350 lbs.	
<b>BK-M20FX-24</b>	20"	6.8 lbs.	23 - 25	2.5"	400 lbs.	

BACK HEIGHT: 22" Tall		FABRIC TYPE: <i>Comfort-Tek</i>			HARDWARE TYPE: <i>Compass</i>	
MODEL NUMBER	HEIGHT	WEIGHT	WC SIZE	THICKNESS	WT. LIMIT	
<b>BK-MTFC-16</b>	22"	5.6 lbs.	15 - 17	2.5"	250 lbs.	
<b>BK-MTFC-18</b>	22"	5.8 lbs.	17 - 19	2.5"	250 lbs.	
<b>BK-MTFC-20</b>	22"	6.2 lbs.	19 - 21	2.5"	300 lbs.	
<b>BK-MTFC-22</b>	22"	6.5 lbs.	21 - 23	2.5"	350 lbs.	
<b>BK-MTFX-24</b>	22"	7.0 lbs.	23 - 25	2.5"	400 lbs.	

24" Wide Back Available with Two Point Mounting Hardware Only

