

# ACTA-RELIEF



Featuring Boa® Technology



Also available as a custom back. See page 52 for details.



FRONT VIEW



SHELL VIEW



WITH HEADREST MOUNTING BRACKET

AVAILABLE  
FALL 2014

# FEATURES

Designed to effortlessly adjust, the Boa® closure system brings comfort and support to the lower and upper thoracic regions.

With the Boa® reel and cable system, equal and stable pressure distribution is accomplished in an "off the shelf" back design.

The open back frame design gives the ability to off-load the pressure typically seen on spinal deformities with up to 4" of relief through the back frame opening.

Head support pad attaches through the use of a 2.5" offset bracket giving the necessary clearances to accommodate more aggressive deformities seen in Kyphotic patients.

Lateral trunk supports (E0956) are optional and available on page 41.

Head supports (E0955) and hardware are optional and available on page 42.

## WEIGHT CAPACITY

16" Wide	250 lbs
18" Wide	250 lbs
20" Wide	300 lbs

## STANDARD SIZES

	16" Wide	18" Wide	20" Wide
18" Height	18 x 16	18 x 18	18 x 20
20" Height	20 x 16	20 x 18	20 x 20
22" Height	22 x 16	22 x 18	22 x 20



## THE BOA® ADVANTAGE

**MICRO-ADJUSTABLE:** Each click of the Boa® reel offers one millimeter of precision fit adjustment.

**SECURE:** Once the Boa® reel is locked into place, it stays that way.

**CUSTOMIZABLE FIT:** Boa® delivers smooth, even adjustment without pressure points for maximum comfort.

**CONVENIENT:** The Boa® closure system offers an adjustment solution free of messy straps, ratchets, and fasteners.

## THE BOA® REEL



Pull to release tension.




Push to engage the system. Turn right to tighten.

## PRESSURE REDISTRIBUTION PROPERTIES


● SEAT CONTACT GAP

● SEAT CONTACT

### TRADITIONAL BACK

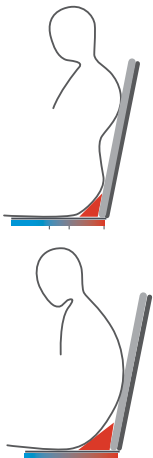

### ACTA-RELIEF BACK

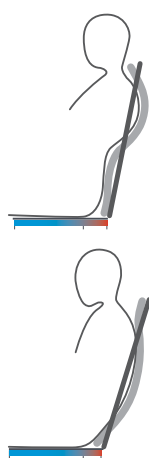

## PART NUMBERS

Hardware	Stretch-Air	Comfort-Tek
COMPASS QUICK RELEASE 2-POINT	BK-KXXSC-YY	BK-KXXFC-YY
COMPASS NON-QUICK RELEASE 2-POINT	BK-KXXSN-YY	BK-KXXFN-YY
COMPASS QUICK RELEASE 4-POINT	BK-KXXSX-YY	BK-KXXFX-YY
COMPASS NON-QUICK RELEASE 4-POINT	BK-KXXSNX-YY	BK-KXXFNX-YY

XX & YY TO BE REPLACED. XX = HEIGHT & YY = WIDTH

## HARDWARE OPTIONS

**C** : Compass Quick Release 2 Point

\*See page 22 for details.


**N** : Compass Non-Quick Release 2 Point

**X** : Compass Quick Release 4 Point


**NX** : Compass Non-Quick Release 4 Point

## COVER OPTIONS

### STRETCH-AIR

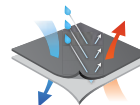

### COMFORT-TEK

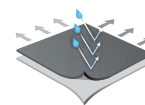

\*See page 3 for details.