


5 Year Warranty

This Non-powered  
Self-adjusting E0373  
Group II Surface is  
Easily Convertible to a  
Powered E0277 APM

NEW  
\*\*\*\*\*1P  
ACK500  
to  
650  
lbs

A

### Heel Slope

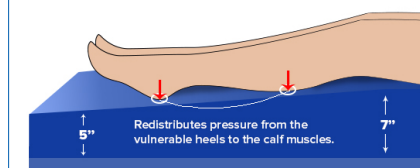

## Balanced Aire Non-Powered Self Adjusting Convertible Mattress

Item #	Product Description	Product Dimensions	Product HCPCS	UOM
BA9600-NP	Non-Powered Self Adjusting Air/Foam Mattress	35" (W) X 80" (L) X 7" (H)	E0373	1/ea
BA9600-P	Powered alternating Pressure Air/Foam Mattress	35" (W) X 80" (L) X 7" (H)	E0277	1/ea
BA9600-NP-84	Non-Powered Self Adjusting Air/Foam Mattress	35" (W) X 84" (L) X 7" (H)	E0373	1/ea
BA9600-P-84	Powered Alternating Pressure Air/Foam Mattress	35" (W) X 84" (L) X 7" (H)	E0277	1/ea
BA9600-NP-42	Non-Powered Self Adjusting Air/Foam Mattress	42" (W) X 80" (L) X 7" (H)	E0373	1/ea
BA9600-P-42	Powered Alternating Pressure Air/Foam Mattress	42" (W) X 80" (L) X 7" (H)	E0277	1/ea

- The (8) 5.5" tri-laminated foam filled air cells are separated into three distinct pressure zones, each zone independently self-adjusting to the patient's needs.
- Designed to prevent and treat any stage pressure ulcer.
- The system is sloped downward in the foot section from 7" to 5" to reduce pressure in the heel zone.
- Add an alternating pressure pump and the Balanced Aire easily converts to a dynamic powered system. (Figure A)
- High Density Foam Perimeter provides solid edge support and facilitates safe patient transfers.

\*Control unit ships only with powered systems.

### specifications

Control Unit Airflow	When Used w/Power Source (P) 4 LPM
Control Unit Cycle Time	When Used w/Power Source (P) 5 Minutes
Control Unit Type	When Used w/Power Source (P) Compressor
Control Unit Alarms	N/A
Mattress Cover Material	Fluid-Resistant Stretch
Weight Capacity	500 lbs.
Warranty	5 Year Limited Mattress   1 Year Control Unit