



110.1 - 150 cm /
43.4 - 59"
18.2 - 50 kg /
40 - 110.2 lbs

Tested in accordance
with the safety standard
FMVSS 213

—
side impact
protection included!

AXKID NEXTKID REHA US






Individual Fit: Perfect adjustment through
flexible side wings and headrest


Thomashilfen



Axkid NEXTKID Reha

- The side wings can be adjusted individually to ensure an optimal seating position.
- A lightweight seat that weighs 8 kg/17.6 lbs for users between 110.1 - 150 cm / 43.4 - 59 inches and 18.2 - 50 kg / 40 - 110.2 lbs.
- The extensive range of accessories not only helps with postural positioning but pairing with the swivel base also helps with ease of transfer in and out of car.
- The height adjustable backrest offers 8 different positions to allow for optimal fit regardless of the child's age or size.

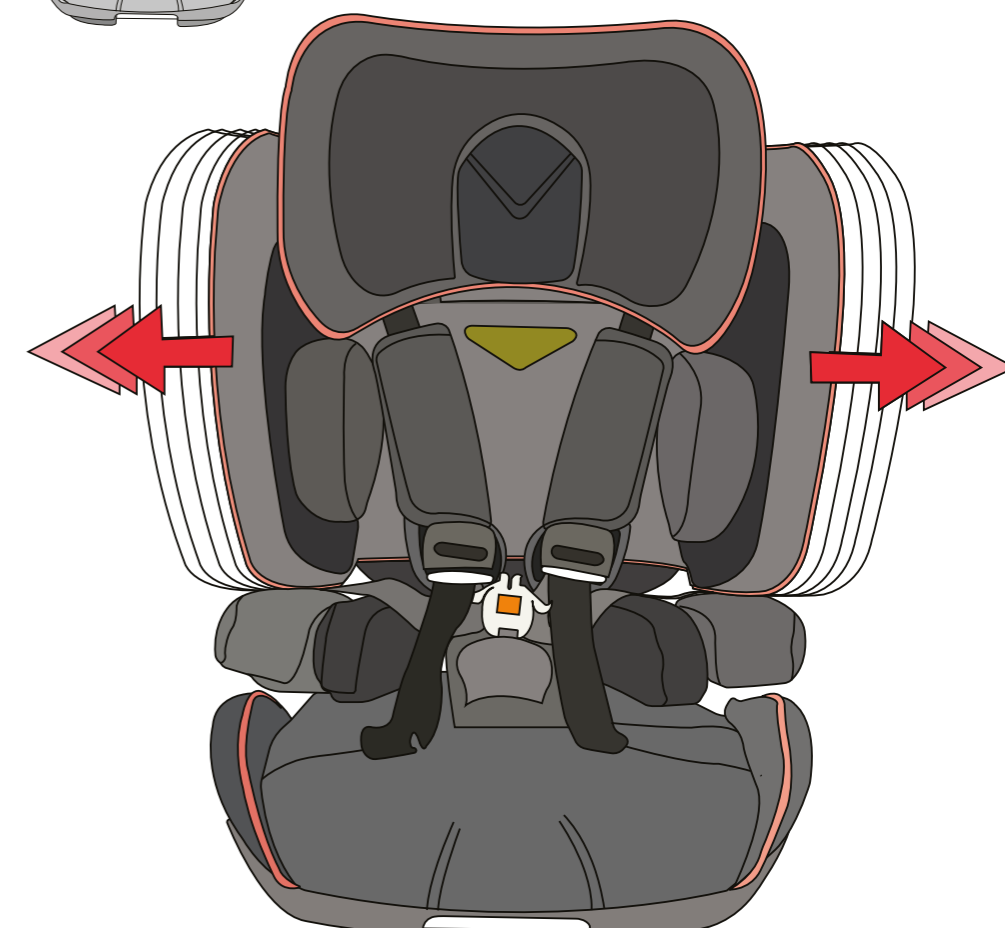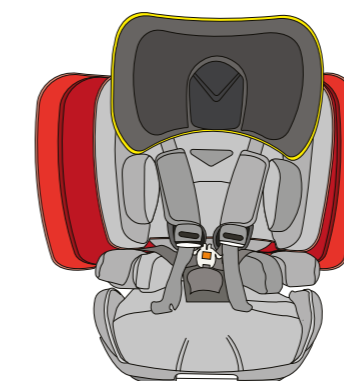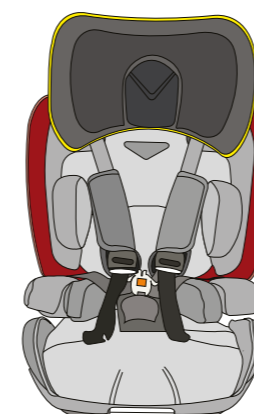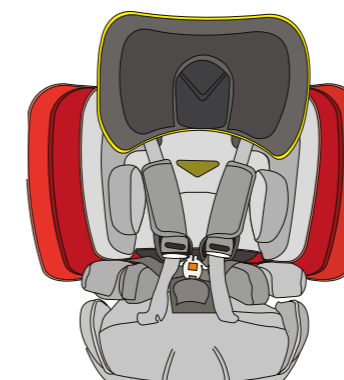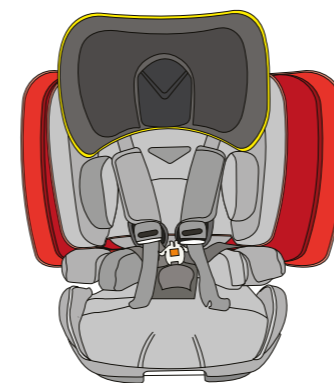


◀ The individually extendable side wings and the separately adjustable headrest adapt perfectly to fit the needs of the user.



Thomashilfen has been supplying rehab car seats for children with special needs for over 35 years. Today, we use our own expertise to develop our own rehab car seats. Incorporating many years of experience in meeting the need for children and their parents.

Gunnar Thomas

Gunnar Thomas, owner

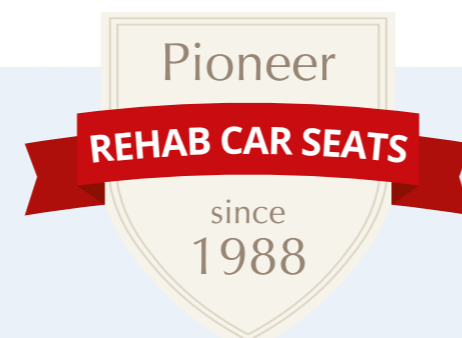

Thomashilfen has the world's largest range of rehab car seats and offers the right support for all ages and disabilities.

BASIC REHAB EQUIPMENT INCLUDED



The basic rehab configuration can optionally be equipped with a suitable swivel and tilting base.



The **SWIVEL** and **TILTING BASE** can be supplemented with a wide range of accessories as required. For easy handling, the release lever is located on the side of the swivel base.

GMFCS
Level
III-V

HARNESS GUIDE

optimizes harness routing and ensures optimal positioning of the safety belt on the pelvis and prevents the belt from slipping off the shoulder or disturbing the neck.

CROTCH PAD

reduces the pressure of the harness buckle on the sensitive abdominal area.

PADDING

The integrated ventilation system reduces heat build-up and sweating.

HEADREST

8-position height adjustable — easily and quickly adjustable, grows step by step with the child. The headrest can be adjusted separately from the side wings.

ISIP SIDE IMPACT PROTECTION (Integrated Side Impact Protection)

PRECISE ADJUSTMENT

for improved upper body posture: large velour upholstery enables precise attachment of lateral trunk support.

The **SIDE WINGS** can be adjusted individually in four stages. Suitable for various disabilities and flexible for many different body shapes and sizes.

ISOFIX / LATCH

Quick and easy installation in just a few seconds. Safety indicator in green/red for correct anchoring. The ISOFIX / Latch bar automatically folds out, making installation even easier (used without swivel and tilt base).

Axkid NEXTKID Reha

► Functions (selection)




LATERAL TRUNK SUPPORT improve upper body support and optimally support by attaching precisely to the velour padding.



The **SIDE WINGS** are individually adjustable, adapt to growth, and ensure an optimal seating position.




LATERAL SEAT UPHOLSTERY provides lateral hip support for smaller users, may be removed for larger users.

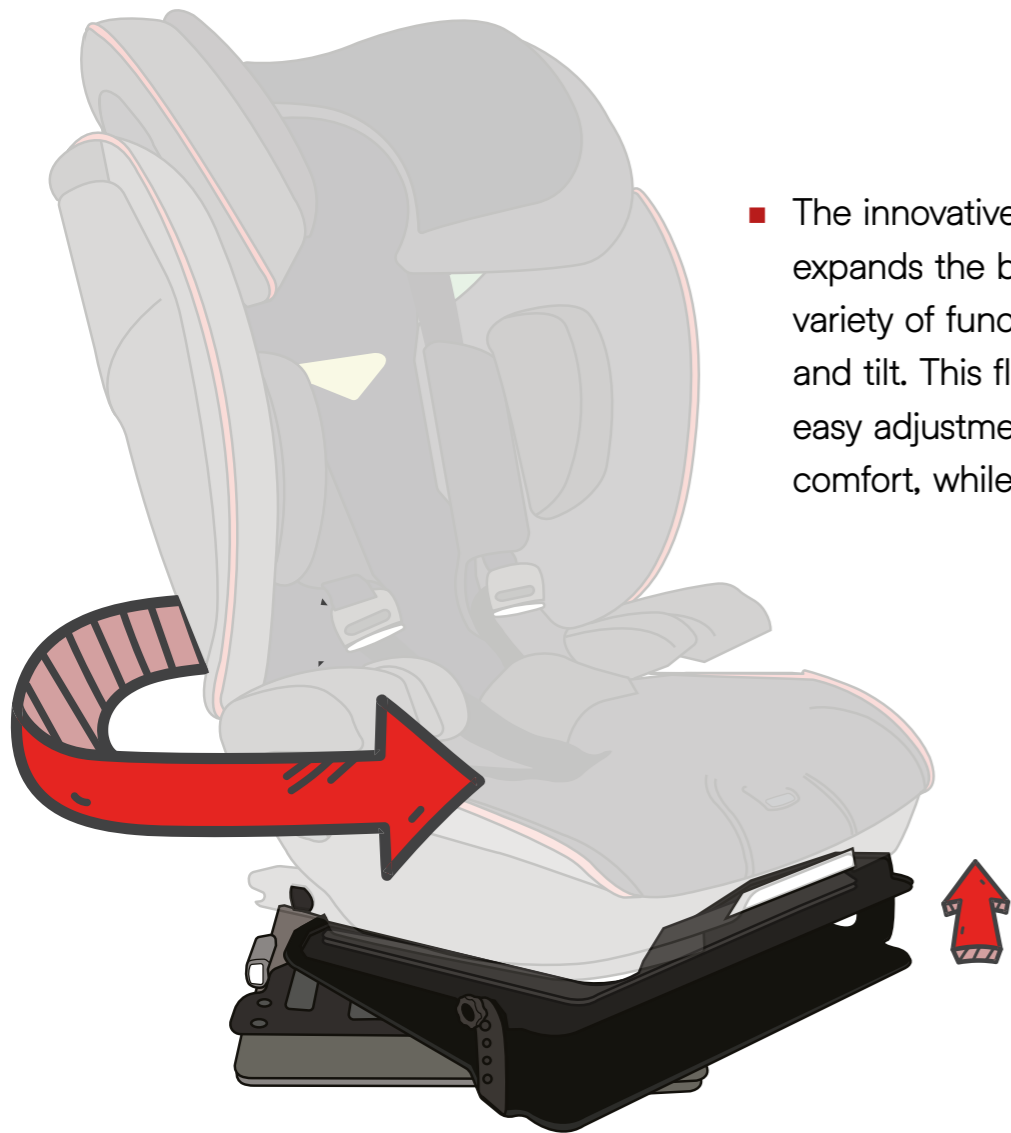


- The innovative high-tech base expands the booster seat with a variety of functions, such as swivel and tilt. This flexibility enables easy adjustment and increases comfort, while keeping the child safe.

HIGH-TECH-BASE COMBINATION

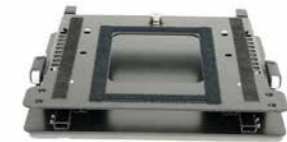

FOOTREST ADAPTER WITH ISOFIX / LATCH

TILTING AND SWIVEL BASE combined for the best possible handling and comfort.

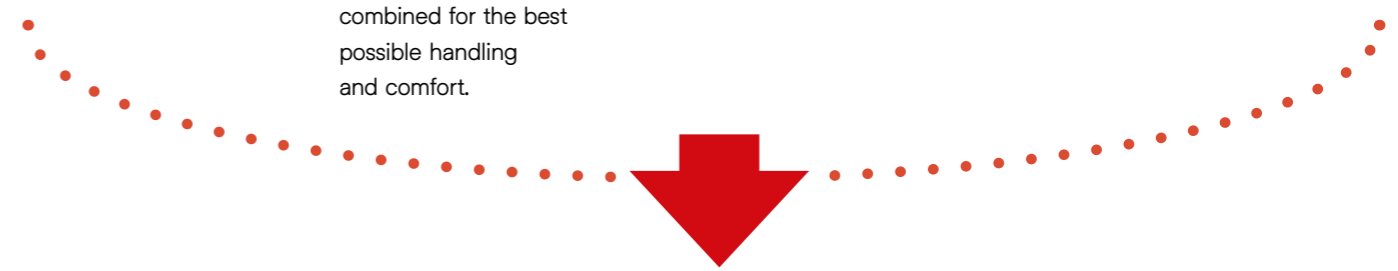

FOOTREST ensures a stable seating position, allowing the child to brace their feet firmly.



SUPPORT LEG ensures stable support of the seat in the car.



ABDUCTION BLOCK (adjustable) improves leg control and prevents the legs from crossing.

► Accessories (selection)



TABLE, ANTERIOR SUPPORT
promotes the child's upright and stable posture, gives the child all-round support and security.



Seat wedge (below)
allows a +10° seat tilt.
To be used without swivel and tilt base.



SEAT DEPTH EXTENSION
offers more leg support and individual adjustment for tall children.

A complete overview of accessories can be found on the order form on the internet: www.thomashilfen.us

Dimensions & weights **Axkid NEXTKID Reha US**

Seat depth (behind crotch strap)	23 cm / 9.1"
Seat depth	37 cm / 14.5"
Seat width	22 - 29 cm / 8.7 - 11.4"
Back height (seat inside)	55 - 75 cm - 21.7 - 29.5"
Shoulder height (harness guide)	32 - 52 cm - 12.6 - 20.4"
Shoulder width	36 - 42 cm / 14.2 - 16.5"
Overall dimensions (w x h x d):	44 x 62 x 48 cm / 17.3 x 24.4 x 18.9"
Weight of seat	8.5 kg / 18.7 lbs
User height (min. - max.)	110.1 - 150 cm / 43.4 - 59"
User weight (min. - max.)	18.2 - 50 kg / 40 - 110.2 lbs

Standard features:

Seat and backrest with integrated 8-position height-adjustable headrest; integrated 5-point belt with central adjuster; height-adjustable shoulder pads; small crotch pad; removable trunk and seat side pads; 4-stage adjustable side wings, adapterplate, chest clip, incl. integrated ISOFIX / Latch connectors, removable covers.


Watch the video:

How quickly and easily the Axkid NEXTKID Reha US is installed in the car.




Rehab car seats from Thomashilfen: The right solution for every journey!




Defender Reha US
(min. - max.)
68.6 - 144.8 cm
27 - 57"
 12 - 54 kg
26.5 - 120 lbs



Hercules Small US
(min. - max.)
110.1 - 150 cm
43.4" - 59"
 18.2 - 50 kg
40 - 110 lbs



Hercules Prime US
(min. - max.)
130 - 180 cm
51.1 - 70.9"
 18.1 - 80 kg
40 - 176.4 lbs

You can find the entire range of rehab car seats at www.thomashilfen.us


Thomashilfen

Thomashilfen North America (ExoMotion, LLC) • 14221 29th St E, #101 • Sumner, WA 98390 • USA
Phone: 866.870.2122 (toll free) • Fax: 866.870.0801 (toll free) • info@thomashilfen.com • www.thomashilfen.us