

Easy to use. Safe for you and your patients.



CareAssist® ES Bed

Hill-Rom

Enhancing Outcomes for Patients and Their Caregivers..

Easy for you

Easy to use, all around the bed.


The CareAssist® ES bed makes every bedside task easier, faster and more efficient.



Point-of-Care® controls
Large, easy-to-read controls just where you want them.



One-button Dining Chair® feature
A single touch of a button moves patient gently in and out of chair position.

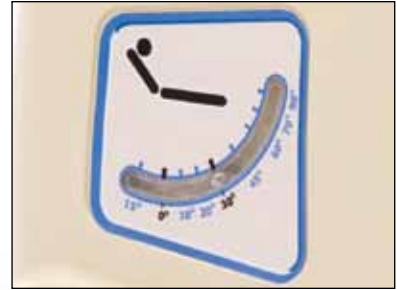

Easy-to-read Line of Site® angle
Monitor patient head angle in any siderail position at point of care.

Built-in bed extender
One caregiver can easily add length to the bed without calling maintenance for an extension.



“Give me a no-hassle bed that saves me time.”



Easy to use, all around the facility.



Built-in connectivity

One of our new smart beds—the CareAssist® ES bed easily connects to your network documenting patient status, saving time and paperwork.



Use almost any surface

Maximum flexibility to care for patients at every level of risk for skin breakdown.



New **Easy to transport**
Now with IntelliDrive® Power Transport, you can transport patients without asking for help.



Battery backup
No need to find an outlet when repositioning the bed during patient transfer.



Fifth wheel
Steer the bed easier with added maneuverability.

An easy choice.

Professional support when you need it

Only Hill-Rom offers access to top-flight clinical professionals and service representatives, who can help you improve protocols for more sustained, measurable improvements in safety.




Safe for you

New Safer transport to protect your back.

The CareAssist® ES bed works hard to create safer patient transport. Its stationary headboard and ergonomically designed handles—at just the right height—make it safer to move patients. And now, with our optional IntelliDrive® feature, patient transport is safer than ever.



“I can move my patient without strain on my back and legs.”

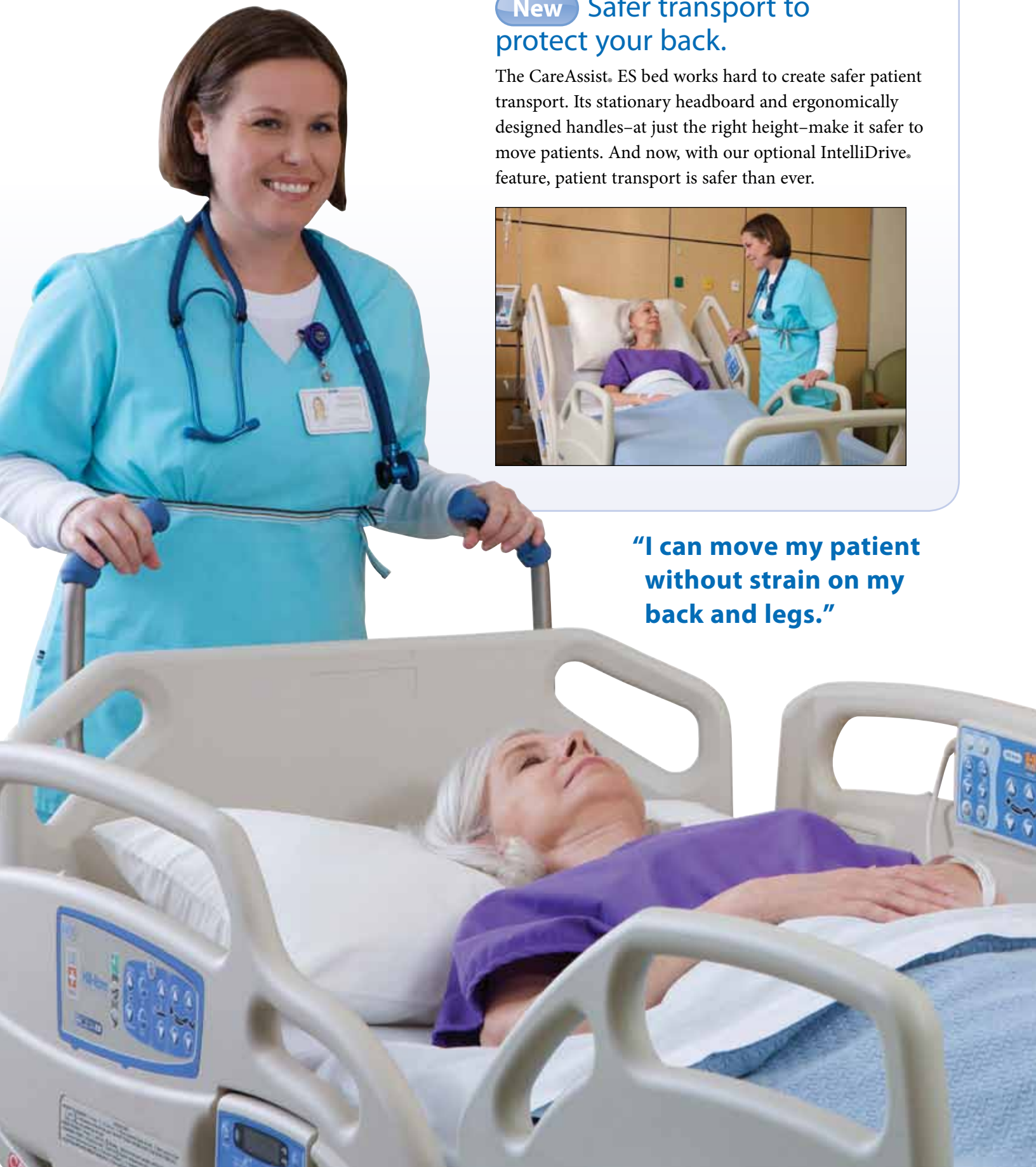

Safer tasks with one-step features.



New Safer handling

Grip the footboard with confidence—it stays firmly in place once it is locked.



Safer positioning

With our electrical vascular foot feature, there's no need to strain to lift the patient's legs and bed frame.



Safer chair positioning

Safely move your patients in and out of chair position, without additional strain on your back.



Safer weighing

Easily weigh your patient in place with superior accuracy—no need to transport to a scale.




Safe for your patients

Protect your patients from skin breakdown.



Safe Skin® program

Led by experienced clinical professionals, our comprehensive program for the prevention of skin breakdown helps improve outcomes that are repeatable and sustainable.

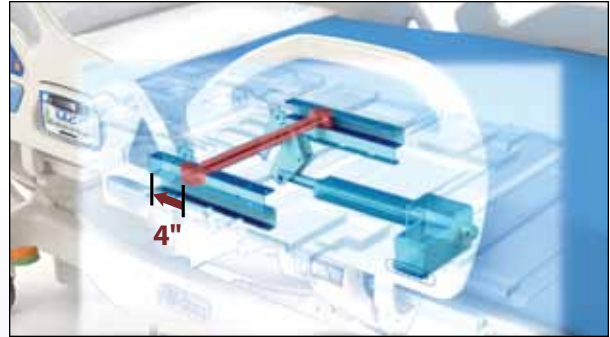

SlideGuard™ Frame Technology


This advanced frame design helps reduce patients from sliding down in bed.

Skin safety comes with your choice of surfaces.



P310 Wound Surface

Multi-zoned prevention surface for moderate acuity patients

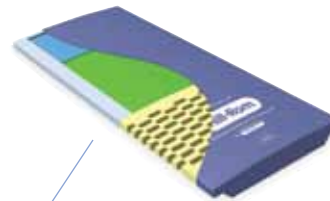

NP200 Wound Surface

Prevention surface with innovative Stimulite® Honeycomb material




Tempur-Pedic® Mattress

The best known name in comfort, now in the medical setting



AccuMax Quantum™ VPC Surface

Non-powered prevention surface with Controlled Release™ technology

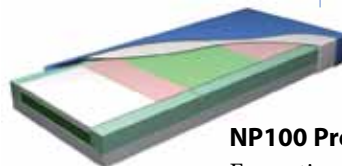

NP100 Prevention Surface

For patients at low risk for skin breakdown

Protect your patients from falls.

Our comprehensive No Falls® program offers evidence-based, continuous improvement of your falls protocols.



SafeView. alerts

Both audible and visible from all angles, these advanced alerts cover bed exit, siderail down, brakes not set and bed height.



Scale-based bed exit monitoring

Highly accurate, sensitive system based on patient weight for fewer false alarms.



Safer brakes


Dual-locking brake casters with an audible alarm if the brake is not set.

“My work is so much easier when I know my patients are safe.”



CareAssist. ES Bed

Easy for you. Safe for you and your patients.

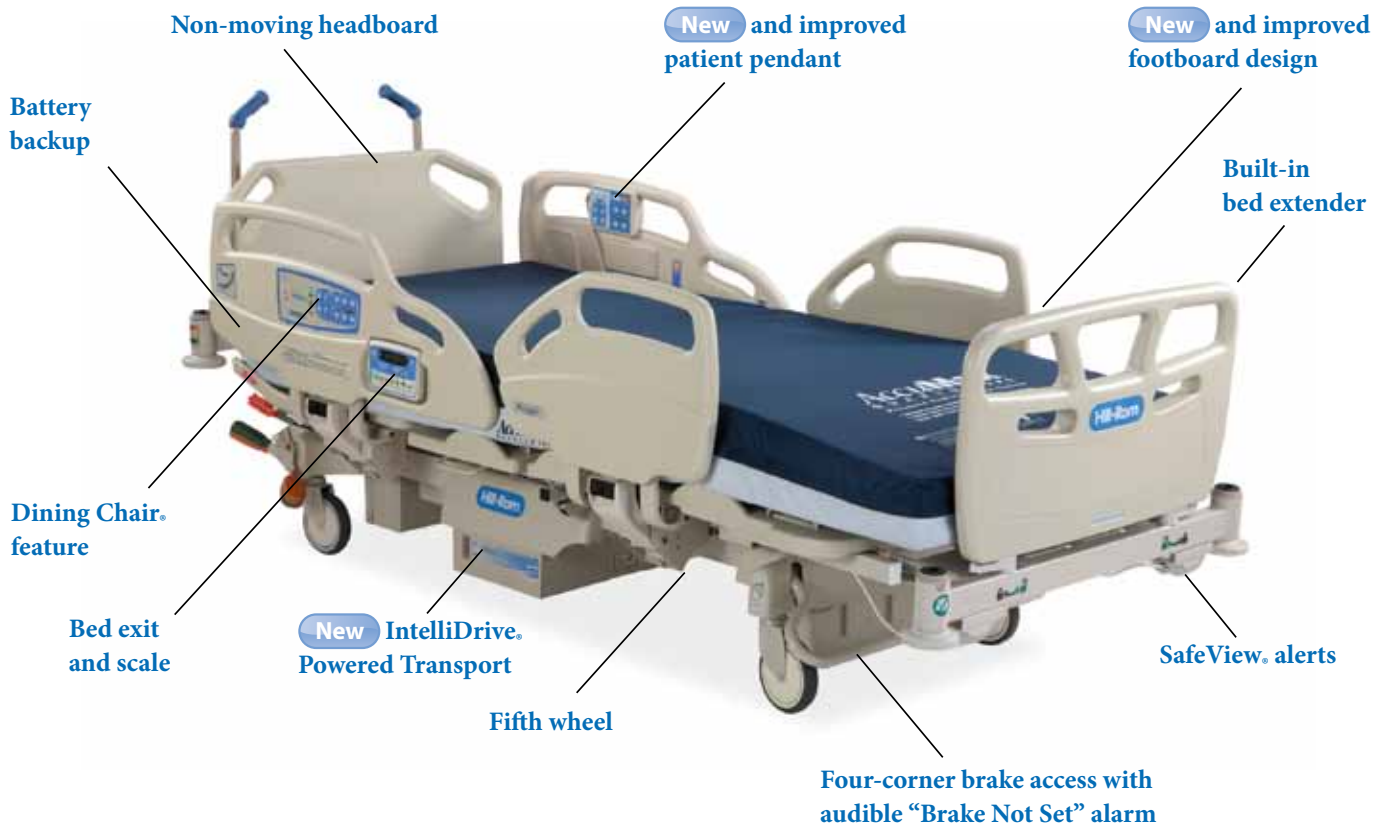

CareAssist. ES Bed Technical Specifications

Hill-Rom. CareAssist ES bed

Frame width	40" (102 cm)
Frame length with head/ footboard (w/ bumpers)	100" (254 cm)
Sleep deck	
width.....	36.00" (91.44 cm)
normal length.....	80.00" (203.20 cm)
extended length	84.00" (213.36 cm)
Bed clearance under frame	5.00" (12.70 cm)

Optional scale system

accuracy	+/- 2.2 lbs or 1.1% of patient weight w/e is greater
weigh in any position	yes
scale repeatability.....	2% of patient weight
Low position floor to deck	15.75" (40 cm)
High position floor to deck	32.5" (82.55 cm)
Max head elevation	65 degrees
Max Trend and reverse Trend	16 degrees
Caster size	6" (15 cm)
Weight limit (Safe working load)	500 lbs

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

© 2011 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.



Enhancing Outcomes for Patients and Their Caregivers.™

USA 800-445-3730
Canada 800-267-2337

www.hill-rom.com

CTG593ra rev 7 6/29/11