

Freestander



Stronger together.

We work with individuals, therapists and carers to design products with both a clinical and an emotional focus. Using the latest research and clinical understanding, we create practical solutions which are easily integrated into family life, because for us, life is about going, enjoying, participating and doing.



24 hour postural care for babies, kids & adults. Sleeping, Sitting, Standing, Walking, Moving, Bathing, Toileting.



Available in 3 sizes, covering a large growth range from 1 - 18 years of age, the Leckey Freestander comes as standard with angle, depth and height adjustable chest, hip and knee supports ensuring the product meets a wide range of clinical needs.

The unique positioning achieved by the pelvic belt enables de-rotation of the hips; alignment of the pelvis and centre of gravity over the base of support; accommodation or correction of flexion at the hips and knees.

The robust design is ideally suited to classroom, therapy room or home environments.

Easy to set up and adjust, making an ideal solution for multiple users. Easy user transfer using standing slings or transfer belts.








Height and angle adjustable activity tray with discreet bowl




A range of head supports is available.




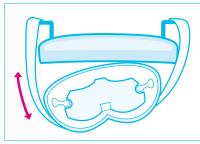

Rotation

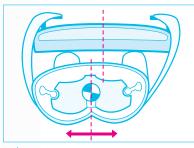

Balance

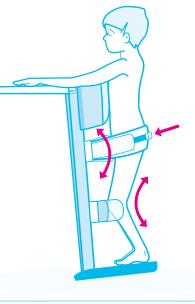

Flexion







The Totstander is a simple and easy to use upright stander. It can be used as an evaluation tool for therapists to assess a child's standing ability.

Promotes active standing, movement and weight shifting depending on the child's level of ability. Available with flexible or rigid chest support, knee supports and sandals which ensure that the product can meet a large range of infants' needs. The Freestander is available in three sizes. The Totstander comes in one size.

Freestander standard set up includes activity tray, bowl and cover, and a choice of: Pelvic belt; knee cups; sandals; flexible laterals. Freestander size 1 and Totstander come as standard in Meadow fabric. Freestander sizes 2 & 3 come in Mondeo Blue.

Product Sizing

Size	1		2		3	
Age	1 - 5		4 - 9 1/2		9 - 18	
Code	LFS/1		LFS/2		LFS/3	
	mm	inches	mm	inches	mm	inches
Freestander						
Chest plate Min Height Max Height	470 680	18.5 27	680 975	27 38	975 1350	38 53
Base Depth Width	650 650	26 26	750 750	30 30	750 750	30 30
Activity Tray Depth Width	400 540	15.7 21.3	425 590	16.7 23.2	470 644	18.5 25.4
User Weight Max	60kg	132lbs	80kg	176lbs	90kg	198lbs
Totstander						
Chest Band Height min max	470 640	18.5 26				
Base Depth Width	450 450	18 18				



Leckey Servicing

Our FREE service and support includes: Free product training Free product set up Free product assessment Free product re-assessment Free repair within 2 year warranty

For service policies on all products outside warranty, please contact Leckey's Customer Service department. James Leckey Design Ltd as manufacturer with sole responsibility declares that all products conform to 93/42/EEC guidelines and EN12182 technical aids for disabled persons general requirements and test methods. Order forms and spare parts lists to extend the service life of the product and allow reissue are available on request or online at www.leckey.com.

Accessories



LFS/#/C Grab rail



LFS/#/L Solid pelvic plate



LFS/#/N Grab posts



1. Contoured head support LIF/2/A,B,C 2. Hyper extension pads LIF/3/A,B 3. Interface frame LIF/1

Suitable for Freestander (not shown) and Pronestander.

Fabric Options



Meadow







Established in 1983, Leckey is a globally recognised pioneer in the research and development of products that help adults and children with disabilities to go, do, enjoy and participate in everyday activities throughout the day and night.

We take a highly clinical approach to product design and development. Through in-depth clinical research studies with leading universities, and extensive trials with occupational therapists, physiotherapists, users and their families, we continue to develop posturally supportive, family friendly products for all day care, at every stage of life.

Through early intervention, childhood and adulthood Leckey's experienced team of designers, therapists and bioengineers work together to develop products that meet the clinical needs of the healthcare professionals and the social needs of the user.

To achieve this, we work with the healthcare professionals, the individuals and carers who use our products everyday. With their help, we create the dependable, durable, proven and high performance products that we are known for worldwide.



















Leckey 19 Ballinderry Road Lisburn BT28 2SA Northern Ireland United Kingdom T: 028 9260 0750 F: 028 9260 0799 E: info@leckey.com W: www.leckey.com



24 hour postural care for babies, kids & adults. Sleeping, Sitting, Standing, Walking, Moving, Bathing, Toileting.