

Jay J2 Deep Contour



SADMERC Code: K0734 adjustable cushion
If other positioning devices (exceeding 1" of postural support) are used a letter from CMS outlining billing practices can be reviewed at www.sunrisemedical.com



POSTURE


- Designed for the client who has moderate to more aggressive postural needs.
- Provides excellent lateral stability as well as forward/rearward stability.
- Provides a stable base of support.
- A firm, anatomically correct deep contoured foam base creates a foundation for postural stability. The foam is closed-cell and extremely lightweight.
- The base and fluid pad can be modified/customized on site to accommodate fixed postures or to correct flexible postures.
 - Modifications can support changing needs over time.
- Optional solid seat inserts provide extra postural stability discouraging the cushion from slinging with the chair upholstery.
- Optional solid drop seat is permanently attached to the bottom of the cushion, replacing wheelchair upholstery, promoting a level-sitting platform when lower seat heights are required.

SKIN

- Most suitable for the client who is at the highest level of risk for skin breakdown or who has very poor skin integrity.
- The 3D Jay Flow™ fluid Tripad™ in combination with the softer foam overlay conforms to bony prominences, distributing pressure more evenly and reducing peak pressures.
- The 3D Jay Flow™ fluid Tripad™ design assists with maintaining fluid under the at-risk bony areas.



Includes an ultra-lightweight, pre-contoured foam base, a pre-contoured deep immersion 3D Jay Flow™ Tripad™, and a zipper enclosed ballistic-edged cover. Optional drop seat can be ordered by placing a "D" at the end of the part number (ex D1208D). To order a Jay J2 with a positioning device installed place a P at the end of the part number."

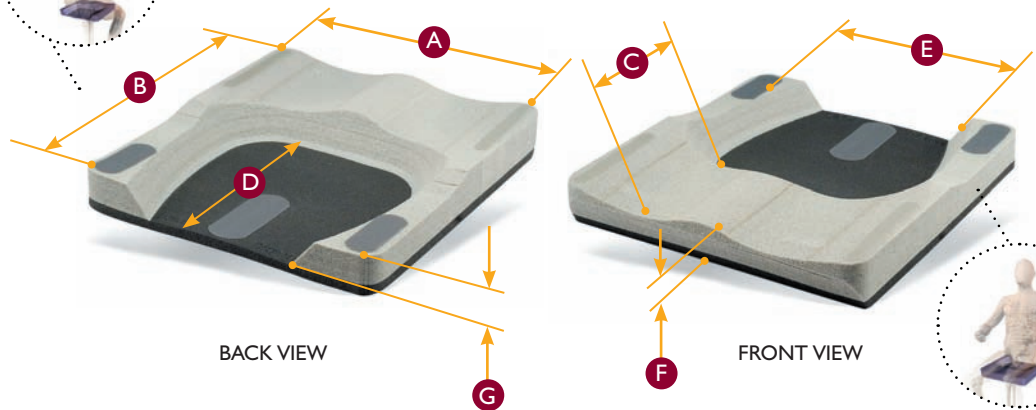

PART #	CHAIR SIZE	ACTUAL WIDTH	ACTUAL LENGTH	FEMORAL SUPPORT LENGTH	PELVIC CONTOUR LENGTH	PELVIC CONTOUR WIDTH	OVERALL HEIGHT	PELVIC CONTOUR DEPTH
		A	B	C	D	E	F	G
2400-MJ	20x20	19.75	19.75	7.25	10.625	18.00	4	3
2406	20x16	19.75	15.75	5.25	8.625	18.50	4	3
2408	20x18	19.75	17.75	6.25	9.750	19.38	4	3
2426	22x16	21.75	15.75	5.25	8.625	20.25	4	3
2428	22x18	21.75	17.75	6.25	9.750	21.63	4	3
2444	14x14	13.75	13.75	4.75	7.625	11.75	4	3
2446	14x16	13.75	15.75	5.75	8.625	12.50	4	3
2448	24x18	23.75	17.75	6.25	9.750	23.50	4	3
2455	15x15	14.75	14.75	5.13	8.000	13.00	4	3
2457	15x17	14.75	16.75	6.00	9.125	13.88	4	3
2460	16x20	15.75	19.75	7.50	10.500	14.25	4	3
2466	16x16	15.75	15.75	5.50	8.500	14.50	4	3
2468	16x18	15.75	17.75	6.38	9.500	15.38	4	3
2475	17x15	16.75	14.75	4.75	8.000	15.25	4	3
2477	17x17	16.75	16.75	5.88	9.125	16.00	4	3
2480	18x20	17.75	19.75	7.50	10.625	16.38	4	3
2486	18x16	17.75	15.75	5.38	8.625	16.50	4	3
2488	18x18	17.75	17.75	6.63	9.500	16.13	4	3

All listed measurements are in inches.