

Jay J2



SADMERC Code: E2613

Includes a back shell, foam back, two lumbar shims, universal mounting hardware and an AirExchange™ cover.

THORACIC

CLINICAL

- ▶ Designed for the client who has minimal to moderate postural needs.
- ▶ Suitable for the client who has changing postural needs over time.
- ▶ Provides excellent posterior pelvis stabilization and posterior thoracic support, facilitating thoracic extension.
- ▶ Addresses the needs of the client with skin integrity issues on the spinous process.

TECHNICAL

- ▶ Lightweight, contoured aluminum shell and contoured foam insert.
- ▶ Angle adjustability of up to 15°.
- ▶ Optional spinal fluid pad insert provides pressure and shear reduction.
- ▶ Height adjustable.
- ▶ Creates minimal loss of seat depth.
- ▶ Shape and size minimizes interference with upper extremity functional activities.

LATERAL

- ▶ Provides mild to moderate lateral support.
- ▶ Optional fixed or swing-away trunk supports ideal for correcting or accommodating a mild to moderate asymmetrical posture.
- ▶ Optional fixed or swing-away trunk supports also ideal for providing symmetrical lateral support.



- ▶ Extended and offset hardware kits are available for a custom fit.
- ▶ Flat and curved pads available.
- ▶ Fluid pads available for the client with skin integrity issues.

Effects to seat depth with relation to recline (standard hooks).

L	M	N
-3°	0.00	1.50
0°	0.42	1.34
3°	0.85	1.18
6°	1.27	1.01
9°	1.70	0.85
12°	2.13	0.69

Effects to seat depth with relation to recline (long hooks).

L	M	N
7°	2.31	0.63
10°	2.74	0.45
13°	3.16	0.27
16°	3.59	0.10
19°	4.01	-0.08
22°	4.44	-0.25

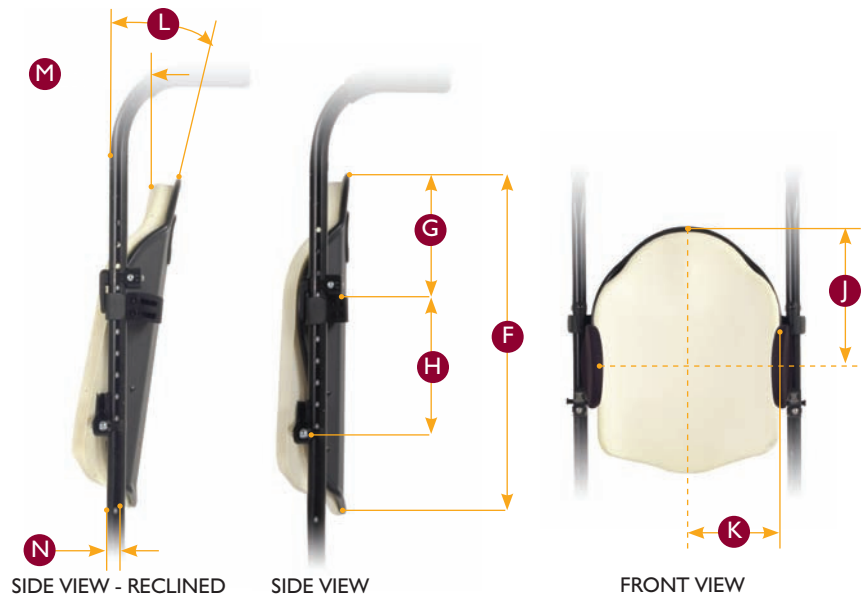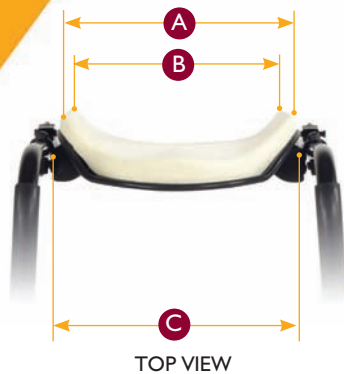

PART #	OUTSIDE WIDTH OF SHELL	WIDTH OF SUPPORT SURFACE	ATTACHING WIDTH OF CHAIRS	UPPER CONTOUR DEPTH OF BACK	LOWER CONTOUR DEPTH OF BACK	HEIGHT OF BACK SHELL	UPPER ATTACHMENT POINT	LOWER ATTACHMENT POINT	VERTICAL LATERAL SUPPORT POSITION	HORIZONTAL LATERAL SUPPORT POSITION
	A	B	C	D	E	F	G	H	J	K
2512K	8.63	8.50	9.25-10.25	1.75	1.250	12.63	3.625	5.25	5.63-9.63	4.00-5.50
2513K	9.50	10.00	10.18-11.13	1.75	1.000	13.00	4.000	5.25	5.63-9.63	4.50-6.00
2514K	11.00	10.75	11.25-12.25	1.75	1.250	14.50	5.125	5.25	5.63-11.63	4.75-6.25
2515K	11.75	12.00	12.25-13.25	2.00	2.125	16.50	6.000	6.00	5.63-13.63	5.25-6.75
2516K	12.88	12.50	14.50-13.25	1.69	2.380	16.75	6.000	6.00	5.50-13.50	5.13-6.63
2517K	14.00	13.50	14.00-15.00	2.25	1.750	16.50	6.000	6.00	5.63-13.63	6.00-7.50
2518K	15.00	14.50	15.25-16.25	2.00	2.750	16.50	6.000	6.00	5.63-13.63	6.50-8.00
2519K	15.75	15.75	16.38-17.38	2.25	2.500	16.63	6.000	6.00	5.63-13.63	7.00-8.50
2520K	17.00	16.00	17.25-18.25	2.25	2.625	16.50	6.000	6.00	5.63-13.63	7.50-9.00

All listed measurements are in inches.