

* HCPCS CODING:
E2607 / E2608

Skin Protection and Positioning



*see chart for dimensions

3D QUADRA GEL - THE QUADRA 3D GEL PACK REDUCES SHEAR FORCE AND HEAT BUILDUP. FOUR 3D COMPARTMENTS ALLOW FOR EVEN DISTRIBUTION OF GEL UNDER THE ISCHIALS AND SEATED AREAS.

COCCYX RELIEF CUTOUT.


TAPERED LEG ADDUCTORS INCREASE IN WIDTH NEAR THE BACK OF THE CUSHION TO GIVE SUPPORT TO THE GREATER TROCHANTER SHELVES.

TAPERED LEG ADDUCTORS, IN CONJUNCTION WITH MEDIAL THIGH SEPARATOR AND DEEP LEG TROUGHS, PROMOTE PROPER LEG ALIGNMENT.

DUAL DENSITY FOAM PROVIDES A SOFT ENTRY LAYER FOR COMFORT AND A MOLDED FOAM BASE LAYER FOR SUPPORT AND PRESSURE DISTRIBUTION.

- STANDARD -

Allows for greater freedom of movement for standard positioning issues.

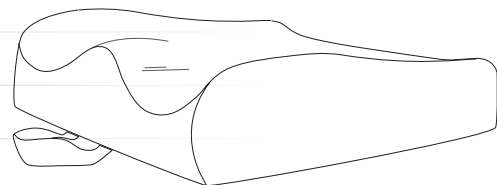

- ANTI-THRUST -

Provides anterior pelvis build-up to prevent forward migration.

- WEDGE -

Controls forward leaning by adding slight anterior height.

-COVERS-

STRETCH-AIR FABRIC — designed for air transmission and heat reduction. This two-ply, breathable laminate combines a high-stretch polyester top layer with a waterproof, breathable vapor barrier.

COMFORT-TEK — fluid proof exterior has supple leathery properties that provide an easily cleaned surface and protection for the inner cushion. Comfort-Tek has multidimensional stretch and helps alleviate pressure on bony prominences.

- STANDARD -

PRODUCT STYLE: <i>Standard</i>		FABRIC TYPE: <i>Stretch-Air</i>		PRODUCT TYPE: <i>Seating Surface</i>		
MODEL NUMBER	W X L	FRONT	REAR	ABDUCTOR	LEG TROUGH	WT. LIMIT
M2-S-1616	16 x 16	3.5"	3"	4.5"	3"	250 lbs.
M2-S-1618	16 x 18	3.5"	3"	4.5"	3"	250 lbs.
M2-S-1620	16 x 20	3.5"	3"	4.5"	3"	250 lbs.
M2-S-1816	18 x 16	3.5"	3"	4.5"	3"	300 lbs.
M2-S-1818	18 x 18	3.5"	3"	4.5"	3"	300 lbs.
M2-S-1820	18 x 20	3.5"	3"	4.5"	3"	300 lbs.
M2-S-2016	20 x 16	3.5"	3"	4.5"	3"	350 lbs.
M2-S-2018	20 x 18	3.5"	3"	4.5"	3"	350 lbs.
M2-S-2020	20 x 20	3.5"	3"	4.5"	3"	350 lbs.
M2-S-2216	22 x 16	3.5"	3"	4.5"	3"	350 lbs.
M2-S-2218	22 x 18	3.5"	3"	4.5"	3"	400 lbs.
M2-S-2220	22 x 20	3.5"	3"	4.5"	3"	400 lbs.
M2-S-2416	24 x 16	3.5"	3"	4.5"	3"	400 lbs.
M2-S-2418	24 x 18	3.5"	3"	4.5"	3"	400 lbs.
M2-S-2420	24 x 20	3.5"	3"	4.5"	3"	400 lbs.

ALL CUSHIONS AVAILABLE WITH RIGID INSERT AND EXTRA INCONTINENT PROOF LINER (PAGE 11)

PRODUCT STYLE: <i>Standard</i>	FABRIC TYPE: <i>Comfort-Tek</i>	PRODUCT TYPE: <i>Seating Surface</i>				
MODEL NUMBER	W x L	FRONT	REAR	ABDUCTOR	LEG TROUGH	WT. LIMIT
M2-F-1616	16 x 16	3.5"	3"	4.5"	3"	250 lbs.
M2-F-1618	16 x 18	3.5"	3"	4.5"	3"	250 lbs.
M2-F-1620	16 x 20	3.5"	3"	4.5"	3"	250 lbs.
M2-F-1816	18 x 16	3.5"	3"	4.5"	3"	300 lbs.
M2-F-1818	18 x 18	3.5"	3"	4.5"	3"	300 lbs.
M2-F-1820	18 x 20	3.5"	3"	4.5"	3"	300 lbs.
M2-F-2016	20 x 16	3.5"	3"	4.5"	3"	350 lbs.
M2-F-2018	20 x 18	3.5"	3"	4.5"	3"	350 lbs.
M2-F-2020	20 x 20	3.5"	3"	4.5"	3"	350 lbs.
M2-F-2216	22 x 16	3.5"	3"	4.5"	3"	350 lbs.
M2-F-2218	22 x 18	3.5"	3"	4.5"	3"	400 lbs.
M2-F-2220	22 x 20	3.5"	3"	4.5"	3"	400 lbs.
M2-F-2416	24 x 16	3.5"	3"	4.5"	3"	400 lbs.
M2-F-2418	24 x 18	3.5"	3"	4.5"	3"	400 lbs.
M2-F-2420	24 x 20	3.5"	3"	4.5"	3"	400 lbs.

ALL CUSHIONS AVAILABLE WITH RIGID INSERT AND EXTRA INCONTINENT PROOF LINER (PAGE 11)

- ANTI-THRUST -

PRODUCT STYLE: <i>Anti-Thrust</i>		FABRIC TYPE: <i>Stretch-Air</i>		PRODUCT TYPE: <i>Seating Surface</i>		
MODEL NUMBER	W X L	FRONT	REAR	ABDUCTOR	LEG TROUGH	WT. LIMIT
M2-SA-1616	16 x 16	4.5"	3.5"	5.5"	4"	250 lbs.
M2-SA-1618	16 x 18	4.5"	3.5"	5.5"	4"	250 lbs.
M2-SA-1620	16 x 20	4.5"	3.5"	5.5"	4"	250 lbs.
M2-SA-1816	18 x 16	4.5"	3.5"	5.5"	4"	300 lbs.
M2-SA-1818	18 x 18	4.5"	3.5"	5.5"	4"	300 lbs.
M2-SA-1820	18 x 20	4.5"	3.5"	5.5"	4"	300 lbs.
M2-SA-2016	20 x 16	4.5"	3.5"	5.5"	4"	350 lbs.
M2-SA-2018	20 x 18	4.5"	3.5"	5.5"	4"	350 lbs.
M2-SA-2020	20 x 20	4.5"	3.5"	5.5"	4"	350 lbs.
M2-SA-2216	22 x 16	4.5"	3.5"	5.5"	4"	350 lbs.
M2-SA-2218	22 x 18	4.5"	3.5"	5.5"	4"	400 lbs.
M2-SA-2220	22 x 20	4.5"	3.5"	5.5"	4"	400 lbs.
M2-SA-2416	24 x 16	4.5"	3.5"	5.5"	4"	400 lbs.
M2-SA-2418	24 x 18	4.5"	3.5"	5.5"	4"	400 lbs.
M2-SA-2420	24 x 20	4.5"	3.5"	5.5"	4"	400 lbs.

ALL CUSHIONS AVAILABLE WITH RIGID INSERT AND EXTRA INCONTINENT PROOF LINER (PAGE 11)

PRODUCT STYLE: <i>Anti-Thrust</i>		FABRIC TYPE: <i>Comfort-Tek</i>			PRODUCT TYPE: <i>Seating Surface</i>	
MODEL NUMBER	W X L	FRONT	REAR	ABDUCTOR	LEG TROUGH	WT. LIMIT
M2-FA-1616	16 x 16	4.5"	3.5"	5.5"	4"	250 lbs.
M2-FA-1618	16 x 18	4.5"	3.5"	5.5"	4"	250 lbs.
M2-FA-1620	16 x 20	4.5"	3.5"	5.5"	4"	250 lbs.
M2-FA-1816	18 x 16	4.5"	3.5"	5.5"	4"	300 lbs.
M2-FA-1818	18 x 18	4.5"	3.5"	5.5"	4"	300 lbs.
M2-FA-1820	18 x 20	4.5"	3.5"	5.5"	4"	300 lbs.
M2-FA-2016	20 x 16	4.5"	3.5"	5.5"	4"	350 lbs.
M2-FA-2018	20 x 18	4.5"	3.5"	5.5"	4"	350 lbs.
M2-FA-2020	20 x 20	4.5"	3.5"	5.5"	4"	350 lbs.
M2-FA-2216	22 x 16	4.5"	3.5"	5.5"	4"	350 lbs.
M2-FA-2218	22 x 18	4.5"	3.5"	5.5"	4"	400 lbs.
M2-FA-2220	22 x 20	4.5"	3.5"	5.5"	4"	400 lbs.
M2-FA-2416	24 x 16	4.5"	3.5"	5.5"	4"	400 lbs.
M2-FA-2418	24 x 18	4.5"	3.5"	5.5"	4"	400 lbs.
M2-FA-2420	24 x 20	4.5"	3.5"	5.5"	4"	400 lbs.

ALL CUSHIONS AVAILABLE WITH RIGID INSERT AND EXTRA INCONTINENT PROOF LINER (PAGE 11)

- WEDGE -

PRODUCT STYLE: <i>Wedge</i>		FABRIC TYPE: <i>Stretch-Air</i>			PRODUCT TYPE: <i>Seating Surface</i>	
MODEL NUMBER	W X L	FRONT	REAR	ABDUCTOR	LEG TROUGH	WT. LIMIT
M2-SW-1616	16 x 16	4.5"	2.5"	5.5"	4"	250 lbs.
M2-SW-1618	16 x 18	4.5"	2.5"	5.5"	4"	250 lbs.
M2-SW-1620	16 x 20	4.5"	2.5"	5.5"	4"	250 lbs.
M2-SW-1816	18 x 16	4.5"	2.5"	5.5"	4"	300 lbs.
M2-SW-1818	18 x 18	4.5"	2.5"	5.5"	4"	300 lbs.
M2-SW-1820	18 x 20	4.5"	2.5"	5.5"	4"	300 lbs.
M2-SW-2016	20 x 16	4.5"	2.5"	5.5"	4"	350 lbs.
M2-SW-2018	20 x 18	4.5"	2.5"	5.5"	4"	350 lbs.
M2-SW-2020	20 x 20	4.5"	2.5"	5.5"	4"	350 lbs.
M2-SW-2216	22 x 16	4.5"	2.5"	5.5"	4"	350 lbs.
M2-SW-2218	22 x 18	4.5"	2.5"	5.5"	4"	400 lbs.
M2-SW-2220	22 x 20	4.5"	2.5"	5.5"	4"	400 lbs.
M2-SW-2416	24 x 16	4.5"	2.5"	5.5"	4"	400 lbs.
M2-SW-2418	24 x 18	4.5"	2.5"	5.5"	4"	400 lbs.
M2-SW-2420	24 x 20	4.5"	2.5"	5.5"	4"	400 lbs.

ALL CUSHIONS AVAILABLE WITH RIGID INSERT AND EXTRA INCONTINENT PROOF LINER (PAGE 11)

PRODUCT STYLE: <i>Wedge</i>	FABRIC TYPE: <i>Comfort-Tek</i>			PRODUCT TYPE: <i>Seating Surface</i>		
MODEL NUMBER	W X L	FRONT	REAR	ABDUCTOR	LEG TROUGH	WT. LIMIT
M2-FW-1616	16 x 16	4.5"	2.5"	5.5"	4"	250 lbs.
M2-FW-1618	16 x 18	4.5"	2.5"	5.5"	4"	250 lbs.
M2-FW-1620	16 x 20	4.5"	2.5"	5.5"	4"	250 lbs.
M2-FW-1816	18 x 16	4.5"	2.5"	5.5"	4"	300 lbs.
M2-FW-1818	18 x 18	4.5"	2.5"	5.5"	4"	300 lbs.
M2-FW-1820	18 x 20	4.5"	2.5"	5.5"	4"	300 lbs.
M2-FW-2016	20 x 16	4.5"	2.5"	5.5"	4"	350 lbs.
M2-FW-2018	20 x 18	4.5"	2.5"	5.5"	4"	350 lbs.
M2-FW-2020	20 x 20	4.5"	2.5"	5.5"	4"	350 lbs.
M2-FW-2216	22 x 16	4.5"	2.5"	5.5"	4"	350 lbs.
M2-FW-2218	22 x 18	4.5"	2.5"	5.5"	4"	400 lbs.
M2-FW-2220	22 x 20	4.5"	2.5"	5.5"	4"	400 lbs.
M2-FW-2416	24 x 16	4.5"	2.5"	5.5"	4"	400 lbs.
M2-FW-2418	24 x 18	4.5"	2.5"	5.5"	4"	400 lbs.
M2-FW-2420	24 x 20	4.5"	2.5"	5.5"	4"	400 lbs.

ALL CUSHIONS AVAILABLE WITH RIGID INSERT AND EXTRA INCONTINENT PROOF LINER (PAGE 11)