

INVACARE®

Matrix® MX1

Back



Yes, you can.

INVACARE®

Matrx® MX1

Back

The latest in back support technology combines aesthetic appeal with superior support and function. The full range of patent-pending hardware options ensures optimal adjustment for the widest range of applications, while the lightweight carbon fiber shell provides exceptional support and durability. Designed for active users who demand peak performance and great looks!



Invacare® Matrx® MX1 Back


Model no. IMX1

Features and Functional Benefits

Carbon Fiber Back Shell:

- Exceptional durability and aesthetic appeal
- Very light weight
- Contoured to provide stability to hips/trunk without restricting upper body mobility

3" CONTOUR DEPTH

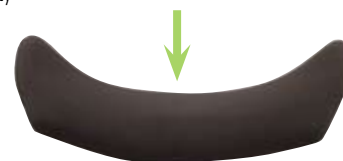

Cover:

- Breathable, attractive spacer mesh fabric
- Cover and foam overlaps back shell for privacy protection and comfort
- Optional privacy flap available



PRIVACY FLAP

Sacral/lumbar support:

- Includes sacral support that is easily adjustable to meet individual needs



SACRAL/LUMBAR SUPPORT

Sizes

Width	14"	15"	16"	17"	18"
Height					
9"	•	•	•	•	•
12"	•	•	•	•	•
16"	•	•	•	•	•
Weight Capacity 250 lb.					

MONO MOUNT



Mounting Hardware Options: (Patent pending)

Mono Mount

- Attaches to wheelchair rigidizer bar providing complete clearance from back canes
- Fully angle adjustable, remains securely fastened in desired position
- Track included provides height adjustment
- Fits rigidizer bars $\frac{3}{4}$ " (19mm) – optional $\frac{7}{8}$ " (22mm) and 1" (25mm) diameter band clamps available
- Allows interface of variety of back model widths

DUAL MOUNT



Dual Mount

- Attaches to wheelchair back canes – requires tools for removal
- Fits back canes 1" (25mm) in diameter – optional $\frac{3}{4}$ " (19mm) and $\frac{7}{8}$ " (22mm) diameter band clamps available
- Attractive, low profile hardware – no sharp edges
- Depth (2"), Angle (30°), Rotation adjustable for optimal positioning
- Width adjustable (2") to fit wheelchairs wider or narrower than back model width, providing optimal fit and support to user *

QUICK RELEASE DUAL MOUNT



Quick Release Dual Mount

- Back cane mount hardware with quick release tabs for easy removal from wheelchair
- Includes band clamps to fit $\frac{3}{4}$ " (19mm), $\frac{7}{8}$ " (22mm) or 1" (25mm) diameter back canes
- Attractive, low profile hardware - no sharp edges
- Depth (2"), Angle (30°), Rotation adjustable for optimal positioning
- Width adjustable (2") to fit wheelchairs wider or narrower than back model width, providing optimal fit and support to user *

*Note – limitation of depth/angle adjustment if used on chairs narrower than model width

Allen McKinley is truly an example of hard work paying off. He is currently employed full time as a Sales Manager, and also works with various companies in the wheelchair industry as a Relationship Manager. When Allen is not working hard in his day jobs, he is spreading smiles at local comedy clubs. He has performed "stand up" for the last three years, pursuing the adrenaline rush that he once enjoyed as a semi-professional snowboarder. Never settling, he has taken his hobby of stand up to 20+ states and performs weekly at the world famous Comedy Store. Understanding that giving time is the greatest gift of all, Allen (with best friend Murphy – also pictured) spends one night every week at Rady's Children's hospital in San Diego running a playroom for patients. He recently received "the Golden Pin" award for his 5+ years of volunteering.

Once Allen understood how much easier great equipment could make his new life in a wheelchair, there was no turning back. Using, the Invacare® Top End® Crossfire™ T7A Ultralightweight Wheelchair with Spinerger FlexRim wheels, Invacare® Matrix® PS Cushion and, of course, the new Invacare® Matrix® MX1 Back, he has the freedom to pursue all kinds of professional avenues to better himself and his life. Sitting in a chair for 12+ hours a day, he knows that to be successful he must invest in the best equipment available.



Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.
Unit 8
Mississauga, Ontario
L4Z 4G4
(800) 668-5324

USA

Motion Concepts
700 Ensminger Dr.,
Suite 112
Tonawanda, New York
14150
(888) 433-6818

Canada

Motion Concepts
84 Citation Dr., Unit 1-7
Concord, Ontario
L4K 3C1
(866) 748-7943
Tel: 905-695-0134
Fax: 905-695-0138

© 2011
Invacare Corporation.
All rights reserved.
Trademarks are
identified by the symbols
™ and ®. All trademarks
are owned by or licensed
to Invacare Corporation
unless otherwise noted.
Matrix is a registered
trademark of Motion
Concepts. Spinerger is a
registered trademark of
Spinerger. Specifications
are subject to change
without notification.

Form No. 11-123
Rev. 09/11

For more information on Invacare® Matrix®
Seating Series Products visit
www.Invacare.com/Rehab

