

Full Platform Exterior Lifts

Base weight lift	92 lbs. ¹
Power chair overall weight	325 lbs.
Maximum product width	28"
Maximum travel distance from the ground	20"
Motor	Sealed 12 volt DC
Overall platform length	49.25" (41" usable surface)
Maximum product wheelbase	35"
Warranty	3-year transferable
Hitch Class	Class II or III ^{2,3}
Vehicle Platform	Full-size van, Mini-van, Truck, Sedan, SUV and Crossover

¹ Weight does not include any accessories

² 190 lbs. maximum product weight with a Class II hitch

³ Vehicle dependent

⁴ Call for details when lifting units with power positioning systems

⁵ Optional lockdown arm required. Call for details.

Full Platform Exterior Lifts

Can easily lift the following current Pride® Power Chairs^{3, 4}

- Jazzy® 600
- Jazzy® 1113 ATS
- Jazzy® Select™
- Jazzy® Select™ 6
- Jazzy® Select™ 6 Ultra
- Jazzy® Select™ 14
- Jazzy® Select™ 14 XL
- Jazzy® Select™ Traveller
- Jet™ 3 Ultra

Can easily lift the following current Pride® Scooters⁵

- Go-Go Elite Traveller® 3 & 4
- Go-Go Elite Traveller® Plus/Plus HD 3 & 4
- Go-Go® Ultra X 3 & 4
- Celebrity® X 3 & 4
- Legend® 3 & 4
- Maxima 3 & 4
- Rally® 3 & 4
- Revo® 3 & 4
- Sundancer® 3
- Victory® 10 3 & 4

Can easily lift the following current Quantum® Power Chairs^{3, 4}

- Quantum® 600/600XL
- Quantum® 610
- Quantum® 1103 Ultra
- Quantum® 1121
- Quantum® 6000
- Quantum® 6000 XL
- Quantum® 6000 Z
- Quantum® R4000

Options

- License plate kit
- Optional battery pack with charger
- Swing-away joint
- Extended swing-away joint
- Power chair weather cover
- Scooter weather cover

¹ Weight does not include any accessories

² 190 lbs. maximum product weight with a Class II hitch

³ Vehicle dependent

⁴ Call for details when lifting units with power positioning systems

⁵ Optional lockdown arm required. Call for details.



Full Platform Exterior Lifts

Your on-the-go lifestyle demands that you bring your power chair with you. Now, Pride Mobility Products® Corp., the number one power chair and scooter manufacturer, offers the lift you need to take your mobility product on the road. Backed by the service you have come to expect and trust from Pride®, with the full support of our expert technical service staff, a Pride® Silver Star™ Lift is the natural complement to your Pride Power Chair.

A Pride Silver Star exterior lift makes it easy to take your Pride Power Chair with you wherever you go, complementing your product perfectly and enhancing your mobility. An all steel frame and protective coating mean Pride Silver Star Lifts are built to be durable, while quality construction means you can count on your lift to be dependable. From loading to lifting, Pride Silver Star Lifts are designed with ease of use in mind.

- 325 lbs. lifting capacity
- 3-year transferable warranty
- Most adaptable exterior lift available
- Drive-on/drive-off platform allows user to enter from either direction
- Lock-down arm available for use with scooters

Notes: The Full Platform Exterior Lift can lift all of the listed power chairs with contoured, Synergy® or Specialty seating. Please call Customer Service when requesting information about lifting a scooter.

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice.

Shock-absorbing
"cradle effect"
protects power
chairs from jarring
vibrations.

New drive-on/
drive-off platform.

Easy to use, three-
point power chair tie
down system allows
for one handed
operation and is
self-retracting and
self-tensioning.

ENHANCED!

All new design dramatically
reduces metal to metal
contact for enhanced
performance and longer
product life.


Innovative, lightweight
self-leveling platform
allows for easy loading
and unloading, even on
inclines or uneven terrain.



Unoccupied lift stowed.



Lift down, ready for product.

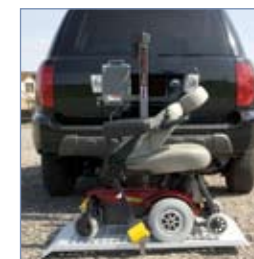

Loaded product secured and
ready for lifting.



Product lifted and ready
for transport.