

INVACARE® TOP END®
Reveal™
Ultralightweight Wheelchair



TOP END®



Yes, you can.®



* Always wear your seat positioning strap. Inasmuch as the seat positioning strap is an option on this wheelchair (you may order with or without the seat positioning strap), Invacare strongly recommends ordering the seat positioning strap as an additional safeguard for the wheelchair user.

INVACARE® TOP END® **Reveal™** Ultralightweight Wheelchair



**Modern design and essential
adjustments with overall
exceptional value.**

* Always wear your seat positioning strap. Inasmuch as the seat positioning strap is an option on this wheelchair (you may order with or without the seat positioning strap), Invacare strongly recommends ordering the seat positioning strap as an additional safeguard for the wheelchair user.

The Invacare® Top End® Reveal™ Wheelchair is the latest ultralightweight performance rigid wheelchair in the Custom Manual family. It combines modern design and essential adjustments with overall exceptional value.

There is a world of possibilities with your Reveal wheelchair. The Reveal Wheelchair is a sleek, super lightweight chair that can easily be adjusted for a customized fit and optimal chair performance. You can tailor the fit to ensure the best wheel access and efficiency by adjusting the center-of-gravity, rear seat height, caster angle, caster placement or back angle. A variety of the most popular options can also be added to the chair to cater to an individual's unique style and needs.

The Reveal wheelchair employs the latest frame material technology, 7005 aluminum, ensuring a smooth and energy-efficient ride. Used in the cycling and aerospace industries, this material helps to increase the chair's overall rigidity and strength without increasing the overall chair weight. This open frame design is super lightweight at a mere 19.71 lb., making the chair more compact than a full frame and easy to transport.

The Reveal wheelchair has it all – the simple adjustments you need, the performance you want and an ultra modern minimalist design – all at an affordable price. Ride in style with the new Invacare® Top End® Reveal™ wheelchair.





**There is
a world of
possibilities
with the
Invacare®
Top End®
Reveal™
wheelchair.**

KEY FEATURES

1 Full offering of tires and wheels including the new standard Lightweight Performance Stylized Spoke Wheel.



2 The stylish grey mesh, no-charge back upholstery offers comfort and a modern look with its soft, breathable material and cushioned inserts around the spine. The grey mesh compliments the new no-charge caster family.




3 The frame allows for ample center-of-gravity and rear seat-to-floor adjustments by using the seat support to fine tune for a more optimal fit.



4 Adopted from Invacare® Top End® popular everyday chairs, the back angle can be adjusted simply by loosening a button-head screw and rotating an adjustment square located on the inside of this back bracket. Adjustments can be made in 4° increments with a total range of 16°.



5 Selection of the most popular Invacare® Matrix® Seating options.

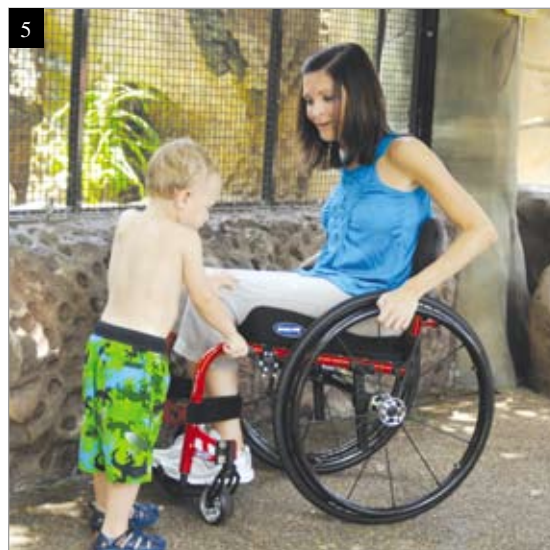

* Always wear your seat positioning strap. Inasmuch as the seat positioning strap is an option on this wheelchair (you may order with or without the seat positioning strap), Invacare strongly recommends ordering the seat positioning strap as an additional safeguard for the wheelchair user.

Specifications



Overall chair width
6.75" plus seat width
with 0° camber



Seat width
12"-20"



Seat depth
12"-20"



Back height
8"-11", 10"-14", 12"-16",
14"-18", 16"-20"



Back angle
98°, 94°, 90°,
86°, 82°



Front frame angle
75°, 85°



Knee-to-heel range
11.25" to 20.25"
(1" increments)*



Center-of-gravity
1", 2", 3"



Front seat-to-floor
15"-20"
(.5" increments)



Rear seat-to-floor
14"-19"
(.5" increments)



Wheel options
22", 24", 25", 26"



Caster options
4", 5", 6"



Weight capacity
250 lb.



Transport weight
13.54 lbs**



Product weight
19.71 lb.***

*Measurement taken from the top of seat rail to bottom of footrest

**Without rear wheels

***16"x16" with 24" lightweight performance spoke wheels, bolt-on seat upholstery, mesh back upholstery with adjustable tension, long seat support and 4" urethane casters.

Frame Colors




Wet Black



Stealth
(Matte) Black



Silver Metallic

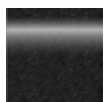

Black Prism



Midnight Blue



Sunny Yellow



Electric Red



Bubblegum
Pink



Grape
Madness



Electric
Blue



Grasshopper
Green



White Snow



Sunkissed
Orange



Tangerine Red

Invacare Corporation

www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

© 2011

Invacare Corporation

All rights reserved.

Trademarks are identified by the symbols TM and [®]. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Matrix is a registered trademark of Motion Concepts. Schwalbe is a registered trademark of RALF BOHLE GmbH Reichshof Germany. Specifications are subject to change without notification.

Form No. 11-146

TOP END[®]



Yes, you can.[®]