

The Invacare® Top End® Pro™ Basketball Wheelchair delivers custom-built performance with adjustments that take only minutes. Perfect for competitors still dialing in their game.



Options

- Clamp-on fixed aluminum side guards
- Spinergy® SLX or Spox wheels
- Sun® Equalizer™ wheels
- 12 tab pushrims
- Kendra tires
- D-Ring positioning belt
- Click straps
- Spoke guards
- Frame protectors
- Alignment gauge
- Double swivel anti-tip bar
- Cushion: soft or hard
- Folding back
- Fold-in side guards

Features

- Custom aluminum frame with offensive 4 bend wing
- Adjustable front and rear seat heights
- Anti-tip/5th wheel: Fixed single swivel adjustable height
- Fixed back
- Adjustable back height
- Adjustable tension back upholstery
- Adjustable Center-of-Gravity
- Fixed camber: 15°, 18° or 20°
- High performance Sun® wheels: 24", 25" or 26"
- Anodized pushrims
- High pressure clincher tires
- Smart fork caster system and super hard 72mm casters
- Platform footrest: Adjustable height, angle and fore/aft
- Detachable and adjustable height swivel anti-tip bar
- Removable offensive wing or bumper
- Powder coat paint



Now Available Per Request!
See a CAD® drawing of your new chair
before it is built.

Top End

Sports and Recreation
Products
4501 63rd Circle North
Pinellas Park, Florida
33781-5914
(727) 522-8677
(800) 532-8677
Fax: (727) 522-1007
www.topendwheelchair.com

Invacare Corporation

www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

©2011

Invacare Corporation.
All rights reserved.
Trademarks are identified
by the symbols ™ and ®.
All trademarks are
owned by or licensed to
Invacare Corporation unless
otherwise noted. Sun is a
registered trademark of
Sun Components Inc.
Spinergy is a registered
trademark of Spinergy, Inc.
AutoCAD is a registered
trademark of Autodesk, Inc.

Form No. 04-188
Rev. 08/11



INVACARE® TOP END®

Basketball Series Wheelchairs



Yes, you can.®

Introducing exciting upgrades in material and design of the Paul Schulte Signature Basketball Chair. The weight was reduced by 2-3 lb. and now it is even faster and more responsive!



New Material: 7000 Aluminum reduces the overall weight of the chair by 2 lb. making it comparable in weight with titanium frames. The big advantage of 7000 aluminum is three-fold. One, the material is approximately 23% stronger than 6061 aluminum. Two, since the material is stronger, we are able to reduce the wall thickness of the aluminum tubes, reducing the overall weight. Best of all, the rigidity and strength of a 7000 aluminum frame is further enhanced by Invacare Top End's sophisticated heat treating and powder coating process.

New Design: The Top End design team has made improvements in multiple components on the Schulte Series BB chairs. The head tube design improvements utilize precision components that maximize stiffness, reduce play and improve overall responsiveness of the casters. Axle receives will now be locked onto the frame using the Nord-Lock bolt securing system. This system minimizes axle play so that a player's rear wheels are secure and solid for a confident feel. Lastly, sleek curved frame members add a sharp look to the camber tube, anti-tips and offensive wing.



Four bend wing. Compact and great for dribbling.



Round wing with round gussets for smooth turning.



Bumper/omit wing for a more defensive setup.




Double amputee front-end option.



Optional adjustable height double swivel anti-tip/fifth wheel.



Optional welded fixed tube sideguards ensure a snug, secure, "one-with-your-chair" feel.

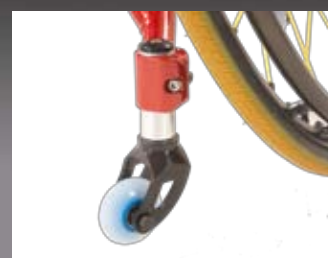

Adjustable height fixed swivel anti-tip with 58mm fifth wheel for the quickest turns ever.




New head tube improvements reduce overall weight by 1 lb. and add precision components that maximize stiffness, reduce play and increase overall performance.

Options

- Adjustable Center-of-Gravity
- Spinergy® SLX or Spox wheels
- Sun® Equalizer™ wheels
- 12 tab pushrims
- Tapered seat width
- Welded seat rigidizer bars w/solid seat
- Double swivel anti-tip/5th wheel
- Click straps
- Double amputee front-end
- Recessed platform footrest
- Adjustable angle footrest
- Clamp-on fixed aluminum side guards
- Fixed side guards w/welded tube
- Fixed aluminum or carbon fiber side guards
- Titanium axles
- Titanium pushrims
- Front frame protectors
- Alignment gauge

Features

- Custom-built fixed frame in 7000 aluminum
- Adjustable height fixed swivel anti-tip with 58mm fifth wheel
- Fixed back w/fixed or adjustable height
- Adjustable tension seat and back upholstery
- Fixed axle position w/fixed wheelbase
- Fixed camber: 15°, 18° or 20°
- High performance Sun® wheels: 24", 25", 26" or 700C
- Anodized pushrims
- High pressure clincher tires
- Precision fork caster system and super hard 72mm casters
- Adjustable height footrest
- Powder coat paint
- Rear frame protector
- Nord lock bolt axle securing system

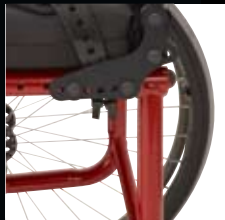

The Paul Schulte signature Chair shown with optional adjustable center-of-gravity.



"When you dream big and give it your all, special things happen."

Paul Schulte



**Custom request? Chair setup questions?
Or just want a computer-generated design
of your chair e-mailed to you before it's
built? No problem!**

Paul Schulte is Top End's design engineer. Along with his degree in Engineering, Paul has achieved numerous championships and MVP's in Juniors, Intercollegiate and the USA men's national team. His years of experience playing with and against the best players in the world are now at your fingertips and he is here to help you up your game.

Contact Paul directly at
pschulte@invacare.com.