



i r i s

intelligent rotation in space



The Quickie IRIS features an innovative design that allows the seat frame to rotate with respect to the base frame, providing a constant center of gravity location while the patient is tilted.



Using a rotational tilt system, Quickie IRIS is unsurpassed in its ability to maintain a user's center-of-gravity in a fixed position. Unlike other systems that incorporate a linear element, the Quickie IRIS' rotational design ensures that both the wheel base and tilt-carriage maintain an almost perfect balance. This design provides the shortest possible wheelbase without compromising stability.



www.sunrisemedical.com




Quickie® IRIS

Advanced— easy to use.

- Unsurpassed maneuverability
- Shortest possible wheelbase without sacrificing stability
- Short, small, light— easy to tilt

Adjustable Carriage

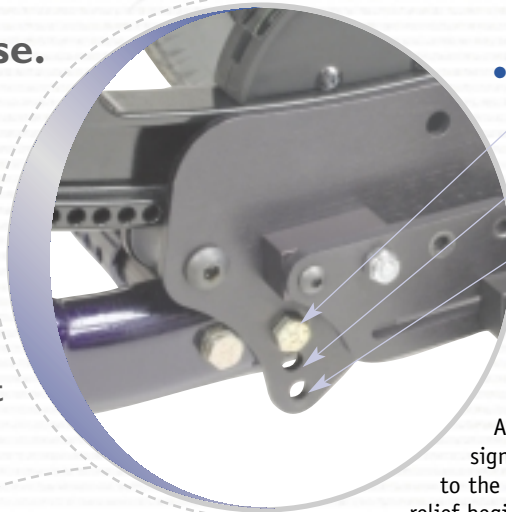

- Tilt Range:
 - 5° to 50°
 - 0° to 55°
 - 5° to 60°

Research studies indicate that greater seat pressure relief is achieved as the seat tilts further backward. Half the load is shifted from the seat to the back at 45°. As the seat tilts further to 60°, a significant portion of the load moves to the backrest. Effective pressure relief begins at 45°.



IRIS SPECIFICATIONS:

| | | |
|-----------------------------------|--|---------------|
| Degrees of Tilt | -5° to 50°, 0° to 55°, 5° to 60° | |
| Seat Width Range | 14" to 22" | |
| Seat Depth Range | 15" to 22" | |
| Seat to floor height | 14.5" to 20.5" | |
| Seat to floor height w/ drop seat | 12.5" to 18.5" | |
| Weight Capacity | Std: | 250 Lbs |
| | Heavy Duty: | 350 Lbs |
| Product Weight* | 35 lbs | |
| Back Height | Angle adj. stroller handle | 16" to 20" |
| | Standard stroller handle | 18", 21", 24" |
| Back Options | Stroller Handle Back Posts (Low, Medium, Tall) | |
| | Height adjustable with stroller handle, & Reclining | |
| Tubing Diameter | Seat rail: | 1" |
| | Back Post: | 1" |
| | Height Adjustable Backrest: | 1" and 7/8" |
| Warranty | Frame: | Lifetime |
| | Accessories: | 1 year |
| Hanger Options | 60°, 70°, 80°, 90° (Swing Away and Heavy Duty), ELR, ALR | |
| Transit Option | Available on Standard and Heavy Duty | |
| Accessories | See Order Form | |
| Overall Length** | 32.5" | |
| Overall Width | Seat width + 8.5" | |

* Without: armrest, footrest, wheels and anti-tip


** Includes anti-tips and footrests

Back Post Bracket



- Depth Adjustment in 1/2" increments
- Easily adjustable— no need to disassemble!
 - Saves time
 - Patent Pending

Modular

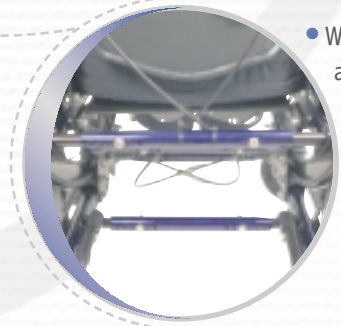

- Width adjustable without adding parts

Adjustable Handle

