



VERSA X



Designed with five adjustable chambers in the rear and a contoured molded base foam in the front. The Versa X is a hybrid air and foam cushion that offers extreme comfort and maximum healing opportunities for any stage of wound.


Liner is made out of a breathable mesh material and allows the full liner with air cells to be washed and dried. Hot or cold settings may be used.

Coccyx, I.T. and sacrum bony prominences each have a dedicated chamber that can be partially or fully off-loaded giving the maximum healing opportunity for any stage of wound.

Independently adjustable lateral chambers in the rear accommodates lateral pelvic asymmetries.

The molded foam cushion is extremely resilient and built with mild positioning features. This design offers the user exceptional stability and comfort while transferring.

The coccygeal bottom insert serves as a pelvic stabilizer and can quickly be adjusted to off-load the high risk, coccygeal and I.T. region through better loading of the greater trochanters.

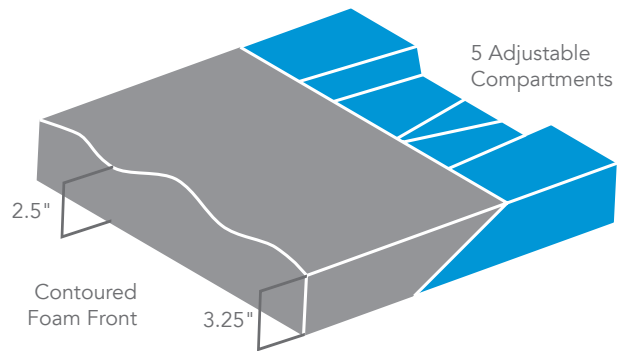

STANDARD SIZES

	16" Wide	18" Wide	20" Wide	22" Wide	24" Wide
16" Deep	16 x 16	18 x 16	20 x 16	22 x 16	24 x 16
18" Deep	16 x 18	18 x 18	20 x 18	22 x 18	24 x 18
20" Deep	16 x 20	18 x 20	20 x 20	22 x 20	24 x 20

WEIGHT CAPACITY

All Sizes	550 lbs
-----------	---------


PART NUMBERS

Cover Style	Stretch-Air	Comfort-Tek
Standard	VRX-S- WWDD	VRX-F- WWDD
With GlideWear®	VRX-SQ- WWDD	VRX-FQ- WWDD


WWDD TO BE REPLACED. WW = WIDTH & DD = DEPTH IN INCHES

COVER OPTIONS

STRETCH-AIR

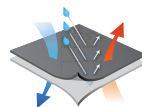

COMFORT-TEK

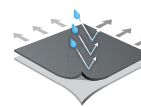

*Available with GlideWear® Technology. See page 3 for details.