

"I sell several [scooter lines], but the scooter that I have used for the past twelve years is a **PaceSaver**."
Bernie Webb, Sec./Tres., Mobility Plus, Inc. (dealer)

Specifications	Eclipse Premier™	Eclipse Atlas™
Overall Length	49.75"	49.75"
Overall Width	24.75"	24.75"
Wheel Base	36.25"	36.25"
Weight		
Rear Drive Unit	41 lbs.	41 lbs.
Base & Tiller	59 lbs.	61 lbs.
Seat Section	24.5 lbs.	29.75 lbs.
Shroud	3.5 lbs.	3.5 lbs.
Maximum Speed	5.40 MPH	5.40 MPH
Range (approximate)	20 - 30 miles	20 - 30 miles
Ground Clearance	4.5"	4.5"
Turning Radius	49"	49"
Pneumatic Tire Size (front/rear)	3" x 10"/3" x 10"	3" x 10"/3" x 10"
Seat Dimensions		
Width	18"	21"
Depth	18"	21"
Back Rest Height	20"	22.5"
Seat Height From Floorpan	16.75" - 18.75"	16.75" - 18.75"
Seat Height With Power Seat	18" - 23"	18" - 23"
Electronic Brakes		
Regenerative, Dynamic	Yes	Yes
Disc Braking	"Posi-Lock" Electric	"Posi-Lock" Electric
Rear Wheel Drive	HD Transaxle	HD Transaxle
Power Source	1 1/4 HP Motor	1 1/4 HP Motor
Structural Capacity	350 lbs	450 lbs
<i>PaceSaver Extra Degree of Protection™</i> - Maximum Incline Rating	12° (22% grade)	12° (22% grade)

Specifications subject to change. Incline stability is determined per US Ansi/Resna Dynamic Stability test 8.2. Stability increases by leaning forward or seat not in highest position.

Five-Year Limited Warranty

Five-year prorated limited transaxle warranty. Three-year prorated limited warranty on Main Frame components. One-year prorated limited warranty covering virtually all internal components, including the electronics.

Lifetime Extended Warranty

Lifetime warranty on frame components. Two-year prorated limited parts warranty. Please consult the Owners Manual for eligibility requirements and details of the warranties. Batteries are not part of either warranty.

Battery Charger - Battery charger's are standard equipment with PaceSavers and are a heavy duty, Micro-Light 24 volt constant current fully automatic charger specially manufactured for Leisure-Lift.

Batteries - PaceSavers are designed to utilize size U-1 sealed deep cycle type batteries. These batteries, supplied by your dealer, may be purchased from Leisure-Lift or supplied by a local source. Leisure-Lift only provides warranty coverage on its own batteries. Because they are 24 volts, PaceSavers require two (2) batteries. Please note: all batteries require some form of maintenance. Deep cycle batteries must be kept charged! Batteries left for some time without using and/or charging will go bad.

Your Professional Local PaceSaver Dealer

PaceSaver® Scooters! Celebrating 21 years of legendary performance.

Color Choices

Ocean Blue

Midnight Black

Candy Apple Red

Velvet Seats

Black

Crimson

Gray

Plush Vinyl Seats

Gray

Tan

Charcoal

Accessories & Options

- U1 Sealed Batteries
- Baskets
- Sliding Seat
- Cane Holder
- Admiral Seat
- Cane Seating
- Review Mirror
- Flat-Free Tires



Flag



Oxygen Holder



Cup Holder



Saddle Bag



Front Bag



Power Seat

Shown with optional Reclining Admiral Seat



Ocean Blue

Please see the Accessories & Options brochure for a complete listing.

PaceSaver®
Scooters

Eclipse™
Premier



- Top Reliability*
- Rock-Solid Stability*
- Beautiful Styling*
- Real Transportability*
- 5 - Year Warranty*
- Maximum Comfort*
- Long-Term Value*

Our workhorse scooter

PaceSaver Scooters!
Celebrating 21 years of legendary performance

Set Your Own Pace!

PaceSaver®

Eclipse™

The Outdoor scooter for a run in the sun.



Bill Burke, PaceSaver

OUR HISTORY

In 1967, Bill and Zella Burke risked all they had to invent and manufacture healthcare products. Mr. Burke said “[I will] dedicate myself to inventing home care products for the infirm”. The resulting family of products allows people to maintain their independence so they can stay in their home close to friends and relatives.

The company that Mr. and Mrs. Burke started employed thirty people; twenty were seniors or disabled and they received national awards for hiring veterans. In our 3rd generation, we continue Mr. Burke’s dedication to building reliable long lasting products and supporting America.

WHAT AMERICAN MADE MEANS TO YOU

Today while many scooters and power chairs are made in low production cost countries like China, PaceSaver continues the 30 year tradition of building to a higher standard of durability and reliability. We are the **only** manufacturer of scooters and power chairs to build all our models in the U.S.A. Because our products are built in the U.S.A., the FDA inspects our plant for design controls, good manufacturing practices and quality systems on a regular basis while foreign plants are not. Industry tests of foreign products reveal poor quality components and designs that do not work as claimed or worse are dangerous if used per the instructions in the owners manual. American made gives you a higher level of Quality Assurance.

AMERICAN MADE MEANS U.S. STOCKED PARTS


U. S. stocked parts keeps your scooter running for years. We stock millions of dollars worth of parts in Kansas City, whether for new product or decades old. Our parts are “over here” not “over seas”! So over 90% of our replacement parts are shipped within one working day.

AMERICAN MADE KEEPS AMERICANS WORKING

Leisure-Lift keeps all of America working! We buy over 5,000 parts manufactured by 580 vendors located in 44 different states. The efforts of people from all over the United States go into our mobility products.



The New Com-For-Back™ Bucket Seat for all day support & comfort.

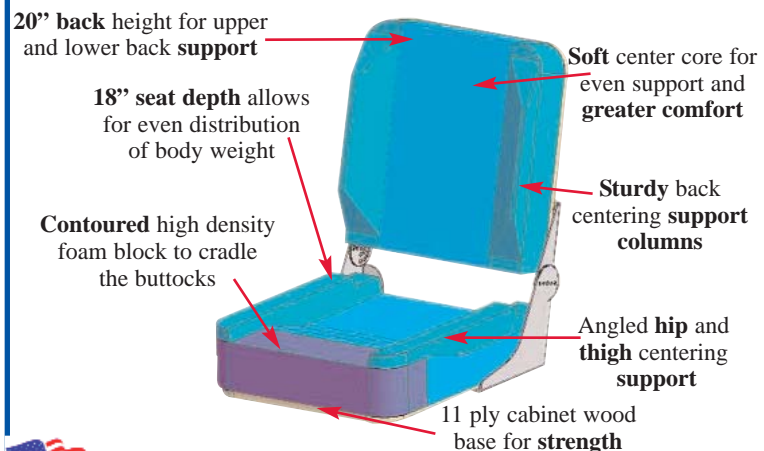

All PaceSaver products are built in America.

Extra Degree of Protection™

- Incline stability is determined per US Ansi/Resna Dynamic Stability Test 8.2.
- A confident ride - indoors & out.

Mirror Glaze™ Lasting Beauty

- Rich “depth of color” reflective finish.
- Color “through and through” ABS substrate.
- Won’t chip, peel or fade.
- “Scratch guard” easy maintenance.
- Deep clear coat hardened acrylic layered.
- Simply wipe with damp rag or wax.
- No paint, No ugly touch ups.

Choices

Our most powerful, rugged workhorse, 4-wheel outdoor/indoor scooter. The stylish, maneuverable Eclipse features **easy access** Soft Touch Delta Tiller for hand and **shoulder comfort**, Com-for-Back high back seat for all day comfort and full support, stylish mag wheels, lights, basket, “Mirror Glaze™” color for **rich long lasting beauty**. The Eclipse has these outstanding performance features: the **industry’s best stability** rating, the **most powerful, easiest to drive**, the **safest state-of-the-art standard controller**, high power, efficient 4 pole motor with heavy duty 3 brake transaxle drive system, a “Smooth ride” front and rear **suspension**, 10” tires for extra outdoor ground clearance, **easy no tool separation** and easy access **on board worldwide charger**. Get out and about on the best American made, most robust 4-wheeler - *Eclipse*.

Eclipse Premier™

Enjoy luxurious full sized comfort and confidence with the Eclipse Premier. Designed for long hill climbs and rugged terrain, Eclipse Premier features 350 lb. weight capacity, 5.40 mph (jogging speed) with 20 to 30 mile range and the industry’s best warranty.

Eclipse Atlas™

The 450 lb. capacity Eclipse Atlas adds a special reinforced frame with 21” wide seat (power seat optional), height, width and length adjustable armrests. Atlas retains the outstanding style, maneuverability, stability, transportability and comfort of the Eclipse.



Shown with optional Admirals Seat

“I sell several [scooter lines], but the scooter that I have used for the past twelve years is a PaceSaver”.

Bernie Webb, Sec./Tres., Mobility Plus, Inc. (dealer)

Digital Electronics

- Powerful 110 amp. controller for more power.
- Push-Too-Fast™ eliminates disengage runaways.
- Power Boost™ technology for inclines.
- Battery Saver™ for extended range.
- Audio-visual indicators provide constant status.



American made controller is enclosed & isolated. VA & CE moisture tested.

- Look under the hood!
- Encased digital controller.
- Isolated from road moisture.
- Double layered moisture protected throttle.

Smooth, Reassuring Control

- Advance Controller Programs for smooth transition through all speed ranges.
- Meets FDA required radio wave runaway protection of 20 volts per meter.

Real Transportability

- Quick & easy single pull-pin separation.
- Separates into 4 easy pieces.
- Fits in most cars, vans and SUV’s.
- Or use our Magnum Power Lift.



Easily Transportable