

CLINICAL

- ▶ Most suitable for the sacral sitter.
- ▶ Designed for the client who has symmetrical posture needing moderate postural support.
- ▶ Provides moderate lateral stability as well as moderate forward/rearward stability.
- ▶ Provides a level base of support while being used on sling upholstery.

TECHNICAL

- ▶ The pre-contoured foam base provides a longer seat well to accommodate fixed “sacral sitting” where the ischials are positioned forward.
- ▶ The reinforced curved bottom counteracts the negative effects associated with sling upholstery providing a level sitting surface.

POSTURE




- ▶ Most suitable for the client at moderate risk of skin breakdown.
- ▶ Ideally the client can make a functional independent weight shift and there is no history of decubitus ulcers.

- ▶ The fluid Tripad™ incorporates left, right and rear sections to maintain the fluid under at-risk bony areas.

SKIN

SADMERC Code: Triad-E2603

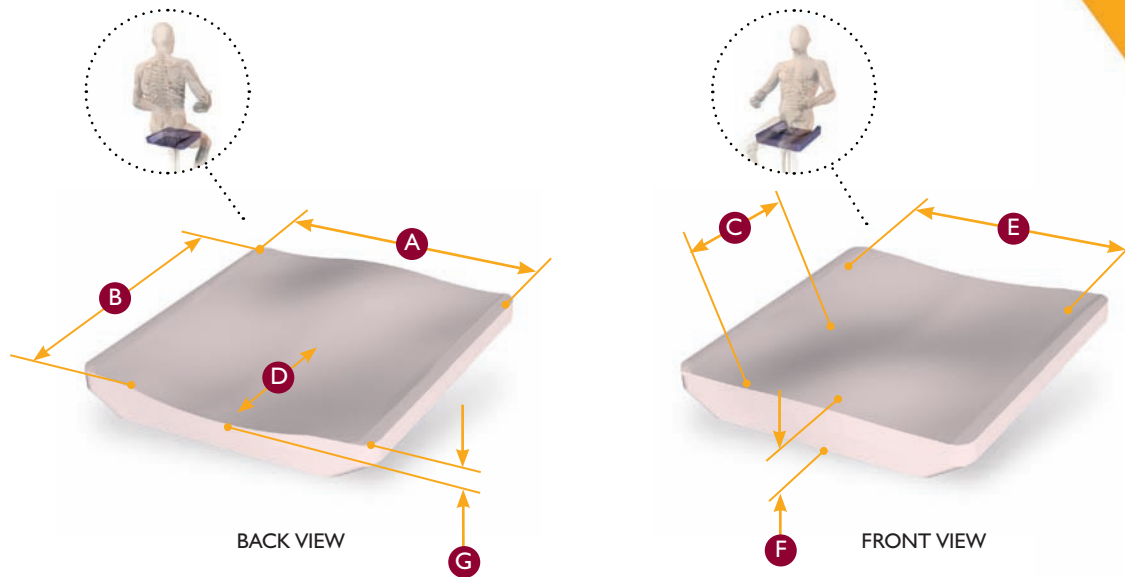

Includes a one-piece cushion, comprised of soft polyurethane foam base, a Jay Flow™ fluid pad for pressure reduction, and a reinforced, no-skid bottom all enclosed by an incontinence-resistant outer casing and covered with a washable stretch cover.

PART #	CHAIR SIZE	ACTUAL WIDTH	ACTUAL LENGTH	FEMORAL SUPPORT LENGTH	PELVIC CONTOUR DEPTH	PELVIC CONTOUR WIDTH	OVERALL HEIGHT	UNLOADED WELL DEPTH
		A	B	C	D	E	F	G
3100B	18x17.5	18.50	17.5	6.00	8.00	17.0	3.75	3.1
3101B	18x19.5	19.25	19.75	7.50	9.00	17.0	3.75	2.9
3105B	16x17.5	17.00	17.75	6.25	8.00	15.0	4.00	2.9
3106B	16x19.5	17.50	19.75	7.00	9.00	15.5	4.00	2.9
3109B	20x19.5	20.75	17.5	6.50	7.75	19.5	3.75	3.1

All listed measurements are in inches.