

## **The Invacare® Top End® Pro™ Tennis Wheelchair ...a totally adjustable Tennis Wheelchair**

For new players, adjustability is key. The **Invacare Top End Pro Tennis Wheelchair** has just that, featuring adjustable front and rear seat heights, center-of-gravity, footrest angle/ height, backrest angle/height and anti-tip. Plus, the **Invacare Top End Pro Tennis Wheelchair** is available in two frame height choices: short or tall. The short frame is ideal for children and the tall frame is for adults, which makes this chair a great choice for programs that have multiple users. Adjusts in minutes, super easy, lightweight and durable. This chair is a winner.



*Invacare® Top End® Pro Tennis Wheelchair shown with optional Spinerger wheels.*

### **Features**

- Adjustable front and rear seat heights
- Anti-tip/5th wheel: Fixed single swivel adjustable height
- Fixed back
- Adjustable back height
- Adjustable tension back upholstery
- Adjustable Center-of-Gravity
- Fixed camber 15°, 18° or 20°
- High performance Sun® wheels: 24", 25", 26"
- Anodized pushrims
- High pressure clincher tires
- Fork caster system and super hard 72mm casters
- Platform footrest: Adjustable height, angle and fore/aft
- Detachable and adjustable height swivel anti-tip bar
- Powder coat paint

### **Options**

- Clamp-on fixed aluminum side guards
- Spinerger® SLX or Spox wheels
- Sun® Equalizer™ wheels
- 12 tab pushrims
- Kenda tires
- Click straps
- Spoke guards
- Frame protectors
- Alignment gauge
- Telescoping anti-tip bar
- Cushion
- Folding Back
- Fold-in side guards

**Invacare Corporation**  
www.invacare.com

**Invacare Top End**  
Sports and Recreation  
Products  
4501 63rd Circle North  
Pinellas Park, Florida  
33781-5914  
(727) 522-8677  
(800) 532-8677  
Fax: (727) 522-1007  
topendwheelchair.com

### **USA**

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### **Canada**

570 Matheson Blvd. E.  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada

(800) 668-5324

© 2011

Invacare Corporation  
All rights reserved.

Trademarks are identified by the symbol ™ and ©. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted.

Spinerger is a registered trademark of Spinerger Inc. Sun Equalizer is a registered trademark of Sun Components. US Open is a registered trademark of United States Tennis Association Inc. AutoCAD is a registered trademark of Autodesk, Inc.

Specifications are subject to change without notification.

Form No. 02-003  
Rev. 08/11



INVACARE® TOP END®

# Tennis Series

Wheelchairs



Yes, you can.®



Now you can get this great chair with exciting upgrades in material and design. We reduced the weight by 2-3 lb. so now it is even faster and more responsive!

**New Material:** 7000 Aluminum reduces the overall weight of the chair by 2 lb. making it comparable in weight with titanium frames. The big advantage of 7000 aluminum is three-fold. One, the material is approximately 23% stronger than 6061 aluminum. Two, since the material is stronger, we are able to reduce the wall thickness of the aluminum tubes, reducing the overall weight. Best of all, the rigidity and strength of a 7000 aluminum frame is further enhanced by Invacare Top End's sophisticated heat treating and powder coating process.

**New Design:** The Top End design team has made improvements in multiple components on the T-5 Tennis Elite chairs. The head tube design improvements utilize precision components that maximize stiffness, reduces play and improves overall responsiveness of the casters. Axle receives will now be locked onto the frame using the Nord-Lock bolt securing system. This system minimizes axle play so that a player's rear wheels are secure and solid for a confident feel. Lastly, sleek curved frame members add a sharp look to the camber tube and anti-tip.




The Invacare® Top End®  
T-5™ Tennis Elite™  
Wheelchair.


## Features

### T-5 Tennis Elite Wheelchair

- Custom-built 7000 aluminum frame
- Fixed swivel anti-tip bar/5th wheel
- Fixed back w/fixed or adjustable height
- Adjustable tension seat and back upholstery
- Fixed axle position w/fixed wheelbase
- Fixed camber 15°, 18° or 20°
- High performance Sun® wheels: 24", 25", 26" or 700C
- Anodized pushrims
- High pressure clincher tires
- 72mm front wheel(s) w/matching mini-fork with inset bearings
- Adjustable height tubular footrest
- Powder coat paint
- Precision fork caster system
- Nord lock bolt axle secure system

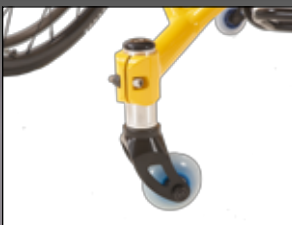


Fixed swivel anti-tip/5th wheel with adjustable height.

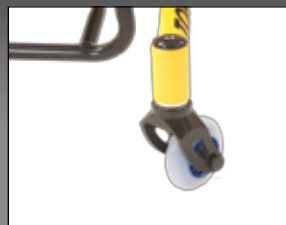


New head tube improvements reduce overall weight by 1 lb and adds precision components that maximizes stiffness, reduces play and increases overall performance



Custom front frame angle, foot space width, axle position and wheelbase.

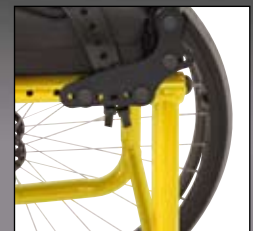

Adjustable center-of-gravity option.

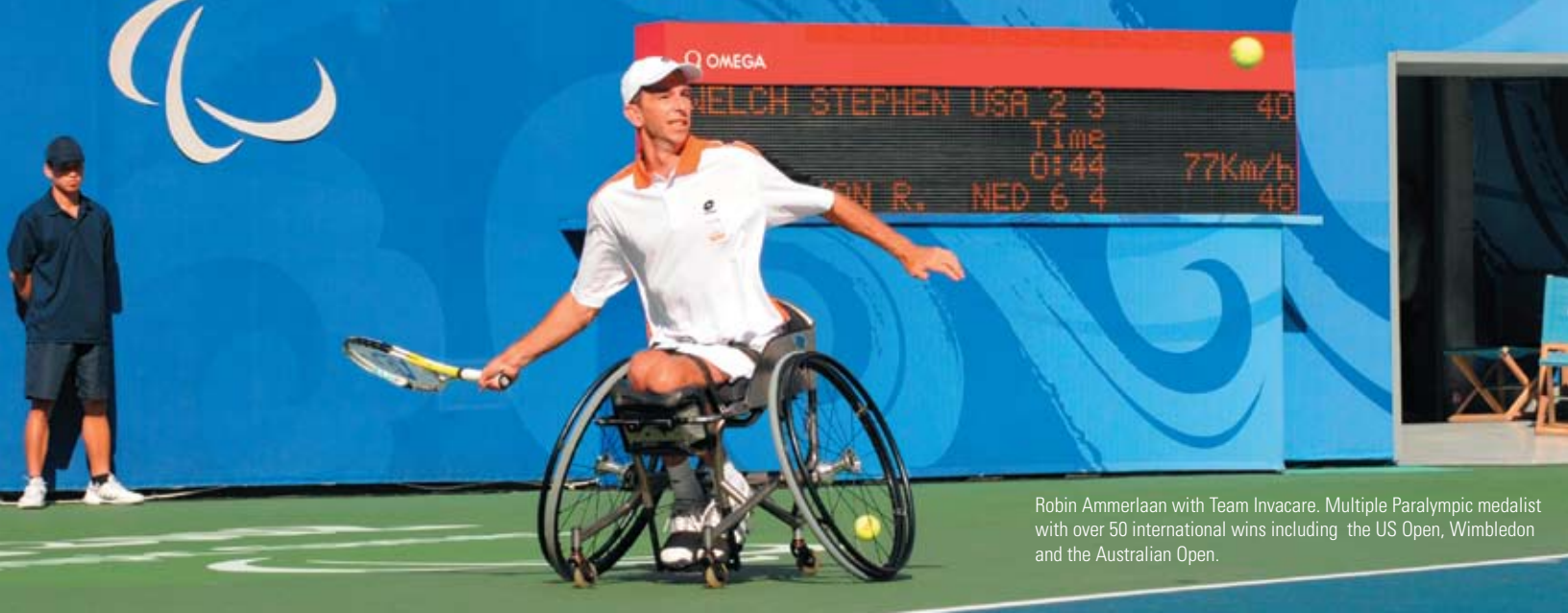

Robin Ammerlaan with Team Invacare. Multiple Paralympic medalist with over 50 international wins including the US Open, Wimbledon and the Australian Open.

## Options

- Spinerger® Spox or SLX wheels
- Sun Fusion wheels
- Adjustable Center-of-Gravity
- Welded seat rigidizer bars w/solid seat
- Telescoping removable anti-tip bar
- Click straps
- Bolt-on fixed side guards
- Fixed side guards w/welded tube
- Fixed aluminum or carbon fiber side guards
- Titanium axles
- Titanium pushrims
- Coated pushrims
- Footrest covers
- Frame protectors
- Custom strapping/toe clips
- Alignment gauge
- Adjustable angle footplate/toe clips
- Recessed platform footplate w/toe stops/footrest strap
- Kenda tires
- Close mount/Cut pushrim tab

The Invacare Top End T-5 Tennis Wheelchair shown with tubular footrest and optional and Spinerger wheels



## Product enhancement continues...

Foot placement is key in tennis. The **T-5 Tennis Elite Wheelchair** is now available with three footrest choices: the standard tubular footrest, adjustable angle footplate which are both adjustable or the platform recessed footrest which has adjustable height, fore/aft and angle with toe stops and footstrap. The T-5 Wheelchair is designed to fit you perfectly and move with minimal effort. To maximize this, many athletes secure themselves with sideguards, click straps and toe clips.

**Now Available Per Request!**  
See a CAD® drawing of your new chair before it is built.

