
Tattoo options

- Mosaic
- 2 Flower garden
- 3 Need for speed



Frame colors

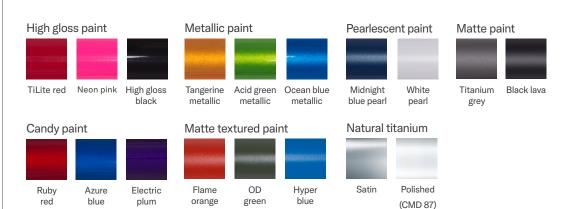

Component colors



















Specifications

TiLite TRA

Frame features	Titanium Dual-Tube Rigid, TiFit, Full adjustability	Front seat width (seat taper)	0", 1", 2", 3"
		Folding seat back	Adj. height 9.5" – 22"
Total chair width*	23"		Adj. angle 80° – 101°
Seat width	12" – 20"	Ergonomic seat	6", 7", 8"
Seat depth	12" – 20"	Transit options	Q-Straint lock down approved
Front seat height	15" – 21"	Transport weight*	11.1 lb
Rear seat height	13" – 21"	Weight limit	265 lb Standard
Front frame angle	70°, 75°, 80°, 85°, or 90°		350 lb Heavy duty
Camber	0°, 2°, 4°, 6°, 8°, or 12°	Warranty	Five years on frame
Center of gravity adjustment	5"		One year on TiLite components

TiLite ZRA

Frame features	Titanium Mono-Tube Rigid, TiFit, Full adjustability
Total chair width*	23"
Seat width	10" – 20"
Seat depth	10" – 20"
Front seat height	16" – 21"
Rear seat height	12.5" – 20"
Front frame angle	70°, 75°, 80°, 85°, or 90°
Camber	0°, 2°, 4°, 6°, 8° or 12°
Center of gravity adjustment	6"

Front seat width (seat taper)	0", 1.5", 2.5"
Folding seat back	Adj. height 9.5" – 22" Adj. angle 80° – 101°
Ergonomic seat	6", 7", 8"
Transit options	Transit tie-down
Transport weight*	10.7 lb
Weight limit	265 lb Standard 300 lb Heavy duty
Warranty	Five years on frame One year on TiLite components

^{*} Weights and widths based on 16" x 16" chair with 0° camber, minimum rear wheel spacing and all standard features. Rear wheels are not included in transport weight. For





TiLite TRA

The dual-tube design of the TRA is one of TiLite's most advanced frame configurations. Designed to minimize weight while also maximizing performance. The adjustability of the center of gravity and backrest enables the user to experience a truly custom fit that meets their needs.

> Custom designed and built to your specifications, with available ergo seating



Find your perfect ride balance with

Dual-tube geometry provides better load distribution and

Adjustable caster fork angle to ensure a precise ride

TiLite ZRA

The classic mono-tube design with a completely adjustable frame. The ZRA is the minamalist look refined.

> Custom designed and built to your specifications, ergo seating available

Adjustable center of gravity to find your perfect ride balance

1.25" diameter tubing to maximize performance and minimize visual impact

> Adjustable caster fork angle to ensure a precision ride

Adjustable seating angle to maximize comfort