



INVACARE® ETUDE™ HC BED



Yes, you can.®

The ETUDE HC Bed features innovative Shear Management Technology (SMT) as a standard feature, reducing shear and friction and providing safety. The SMT feature is provided through the auto-regression of the head section while repositioning in ETUDE HC Bed through use of the eight function hand control.

Features and Options

Modern and Flexible

The ETUDE HC Bed delivers a modern, streamlined and discreet design perfectly complementing any interior. The ETUDE HC Bed can be configured with attractive wood side panels and a choice of safety rails to accommodate specific transfer, safety and support requirements. The ETUDE HC Bed provides a deck height ranging from 7.8" to 26.5". The simple to use eight function hand control enables the user or caregiver to lock out specific functions for increased safety.

Easy Handling & Storage

The ETUDE HC Bed recognizes that portability, storage and set-up are important features for users and caregivers alike. The simplistic design allows for simple assembly without tools and transport can be managed safely as each section of ETUDE HC Bed weighs less than 52 lbs. Smooth roll 3" locking casters enable the bed to be easily maneuvered.

Highest Safety levels

The ETUDE HC Bed conforms to the IEC bed standard IEC 60601-2-52. This equates to the highest degree of safety for both the client and the caregiver. The eight function pendant allows the caregiver to lock out any of the bed's functions for enhanced patient safety.

***HCPCS Code: EZ0266 - PDAC Pending**



Width

Outside:
40" / 102 cm
Inside:
36" / 90 cm



Length

Outside:
86" / 219 cm
Inside:
82" / 200 cm



Height adjustable

Low:
7.9" - 23.5" / 20-60 cm
High:
10.6" - 26.4" / 27-67 cm



Mattress Support

33.5" - 9.5" - 14" - 22"



Backrest Angle

0 - 70°



Thigh Angle

0 - 18°



Legrest Lift

0 - 25°



Tilt Option

11°



Product Weight

183 lbs. /
83 kg



Max Weight Capacity






Bed
440 lbs. / 200 kg
User:
360 lbs. / 165 kg

Limited Warranty

5 years on frame, 3 years on actuators, 2 years on pendant and electronics.

*All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further details.

Invacare beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new and unused after August 1, 2007, meet the voluntary FDA Guidelines for bed rail entrapment. Visit www.invacare.com/bedsafety for additional safety and entrapment information.

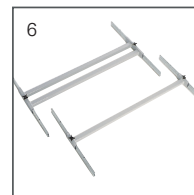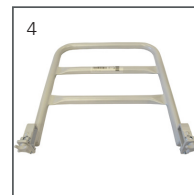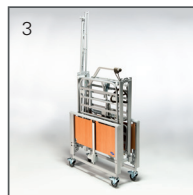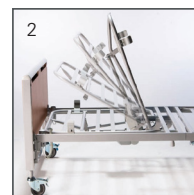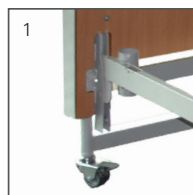

1. Built in Dual Height Brackets
2. Shear Management Technology
3. Transport and Storage Kit
4. Support Handle
5. Lifting Support Trapeze Kit
6. 4" and 8" Extension Kit

Frame and Color of wood work

Aluminum Gray



Beech



Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4 Canada
(800) 668-5324

© 2019

Invacare Corporation.

All rights reserved.

Trademarks are identified by the symbols TM and [®]. All trademarks are owned by or licensed to Invacare Corporation or its subsidiaries unless otherwise noted.

Specifications are subject to change without notification.

Form No. 19-115

Rev. 10/19

191015



Yes, you can.®